

Factor IX19. GlycoPEGylation of Factor IX produced in CHO cells

This example sets forth the preparation of asialoFactor IX and its sialylation with CMP-sialic acid-PEG.

- 5 **Desialylation of rFactor IX.** A recombinant form of Coagulation Factor IX (rFactor IX) was made in CHO cells. 6000 IU of rFactor IX were dissolved in a total of 12 mL USP H₂O. This solution was transferred to a Centricon Plus 20, PL-10 centrifugal filter with another 6 mL USP H₂O. The solution was concentrated to 2 mL and then diluted with 15 mL 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, 5 mM CaCl₂, 0.05% NaN₃ and then reconcentrated.
- 10 The dilution/concentration was repeated 4 times to effectively change the buffer to a final volume of 3.0 mL. Of this solution, 2.9 mL (about 29 mg of rFactor IX) was transferred to a small plastic tube and to it was added 530 mU α 2-3,6,8-Neuraminidase-agarose conjugate (*Vibrio cholerae*, Calbiochem, 450 μ L). The reaction mixture was rotated gently for 26.5 hours at 32 °C. The mixture was centrifuged 2 minutes at 10,000 rpm and the supernatant
- 15 was collected. The agarose beads (containing neuraminidase) were washed 6 times with 0.5 mL 50 mM Tris-HCl pH 7.12, 1 M NaCl, 0.05% NaN₃. The pooled washings and supernatants were centrifuged again for 2 minutes at 10,000 rpm to remove any residual agarose resin. The pooled, desialylated protein solution was diluted to 19 mL with the same buffer and concentrated down to ~2 mL in a Centricon Plus 20 PL-10 centrifugal filter. The
- 20 solution was twice diluted with 15 mL of 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, 0.05% NaN₃ and reconcentrated to 2 mL. The final desialylated rFactor IX solution was diluted to 3 mL final volume (~10 mg/mL) with the Tris Buffer. Native and desialylated rFactor IX samples were analyzed by IEF-Electrophoresis. Isoelectric Focusing Gels (pH 3-7) were run using 1.5 μ L (15 μ g) samples first diluted with 10 μ L Tris buffer and mixed with 12 μ L
- 25 sample loading buffer. Gels were loaded, run and fixed using standard procedures. Gels were stained with Colloidal Blue Stain (Figure 154), showing a band for desialylated Factor IX.

- Preparation of PEG (1 kDa and 10 kDa)-SA-Factor IX.** Desialylated rFactor-IX (29 mg, 3 mL) was divided into two 1.5 mL (14.5 mg) samples in two 15 mL centrifuge
- 30 tubes. Each solution was diluted with 12.67 mL 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, 0.05% NaN₃ and either CMP-SA-PEG-1k or 10k (7.25 μ mol) was added. The tubes were

inverted gently to mix and 2.9 U ST3Gal3 (326 μ L) was added (total volume 14.5 mL). The tubes were inverted again and rotated gently for 65 hours at 32 °C. The reactions were stopped by freezing at -20 °C. 10 μ g samples of the reactions were analyzed by SDS-PAGE. The PEGylated proteins were purified on a Toso Haas Biosep G3000SW (21.5 x 30 cm, 13 μ m) HPLC column with Dulbecco's Phosphate Buffered Saline, pH 7.1 (Gibco), 6 mL/min. The reaction and purification were monitored using SDS Page and IEF gels. Novex Tris-Glycine 4-20% 1 mm gels were loaded with 10 μ L (10 μ g) of samples after dilution with 2 μ L of 50 mM Tris-HCl, pH 7.4, 150 mM NaCl, 0.05% NaN₃ buffer and mixing with 12 μ L sample loading buffer and 1 μ L 0.5 M DTT and heated for 6 minutes at 85 °C. Gels were stained with Colloidal Blue Stain (Figure 155) showing a band for PEG (1 kDa and 10 kDa)-SA-Factor IX.

20. Direct Sialyl-GlycoPEGylation of Factor IX

This example sets forth the preparation of sialyl-PEGylation of Factor IX without prior sialidase treatment.

Sialyl-PEGylation of Factor-IX with CMP-SA-PEG-(10 kDa). Factor IX (1100 IU), which was expressed in CHO cells and was fully sialylated, was dissolved in 5 mL of 20 mM histidine, 520 mM glycine, 2% sucrose, 0.05% NaN₃ and 0.01% polysorbate 80, pH 5.0. The CMP-SA-PEG-(10 kDa) (27 mg, 2.5 μ mol) was then dissolved in the solution and 1 U of ST3Gal3 was added. The reaction was complete after gently mixing for 28 hours at 32°C. The reaction was analyzed by SDS-PAGE as described by Invitrogen. The product protein was purified on an Amersham Superdex 200 (10 x 300 mm, 13 μ m) HPLC column with phosphate buffered saline, pH 7.0 (PBS), 1 mL/min. R_t = 9.5 min.

Sialyl-PEGylation of Factor-IX with CMP-SA-PEG-(20 kDa). Factor IX (1100 IU), which was expressed in CHO cells and was fully sialylated, was dissolved in 5 mL of 20 mM histidine, 520 mM glycine, 2% sucrose, 0.05% NaN₃ and 0.01% polysorbate 80, pH 5.0. The CMP-SA-PEG-(20 kDa) (50 mg, 2.3 μ mol) was then dissolved in the solution and CST-II was added. The reaction mixture was complete after gently mixing for 42 hours at 32°C. The reaction was analyzed by SDS-PAGE as described by Invitrogen.

The product protein was purified on an Amersham Superdex 200 (10 x 300 mm, 13 μ m) HPLC column with phosphate buffered saline, pH 7.0 (Fisher), 1 mL/min. R_t = 8.6 min.

21. Sialic Acid Capping of GlycoPEGylated Factor IX

This examples sets forth the procedure for sialic acid capping of sialyl-glycoPEGylated peptides. Here, Factor-IX is the exemplary peptide.

Sialic acid capping of N-linked and O-linked Glycans of Factor-IX-SA-PEG (10 kDa). Purified r-Factor-IX-PEG (10 kDa) (2.4 mg) was concentrated in a Centricon® Plus 20 PL-10 (Millipore Corp., Bedford, MA) centrifugal filter and the buffer was changed to 50 mM Tris-HCl pH 7.2, 0.15 M NaCl, 0.05% NaN_3 to a final volume of 1.85 mL. The protein solution was diluted with 372 μ L of the same Tris buffer and 7.4 mg CMP-SA (12 μ mol) was added as a solid. The solution was inverted gently to mix and 0.1 U ST3Gal1 and 0.1 U ST3Gal3 were added. The reaction mixture was rotated gently for 42 hours at 32 °C.

A 10 μ g sample of the reaction was analyzed by SDS-PAGE. Novex Tris-Glycine 4-12% 1 mm gels were performed and stained using Colloidal Blue as described by Invitrogen. Briefly, samples, 10 μ L (10 μ g), were mixed with 12 μ L sample loading buffer and 1 μ L 0.5 M DTT and heated for 6 minutes at 85 °C (Figure 156, lane 4).

Factor VIIa

22. GlycoPEGylation of Recombinant Factor VIIa produced in BHK cells

This example sets forth the PEGylation of recombinant Factor VIIa made in BHK cells.

Preparation of Asialo-Factor VIIa. Recombinant Factor VIIa was produced in BHK cells (baby hamster kidney cells). Factor VIIa (14.2 mg) was dissolved at 1 mg/ml in buffer solution (pH 7.4, 0.05 M Tris, 0.15 M NaCl, 0.001 M CaCl_2 , 0.05% NaN_3) and was incubated with 300 mU/mL sialidase (*Vibrio cholera*)-agarose conjugate for 3 days at 32 °C. To monitor the reaction a small aliquot of the reaction was diluted with the appropriate buffer and an IEF gel performed according to Invitrogen procedures (Figure 157). The mixture was centrifuged at 3,500 rpm and the supernatant was collected. The resin was washed three times (3x2 mL) with the above buffer solution (pH 7.4, 0.05 M Tris, 0.15 M NaCl, 0.05% NaN_3) and the combined washes were concentrated in a Centricon-Plus-20. The remaining

solution was buffer exchanged with 0.05 M Tris (pH 7.4), 0.15 M NaCl, 0.05% NaN₃ to a final volume of 14.4 mL.

Preparation of Factor VIIa-SA-PEG (1 kDa and 10 kDa). The desialylation of Factor VIIa solution was split into two equal 7.2 ml samples. To each sample was added either CMP-SA-5-PEG(1 kDa) (7.4 mg) or CMP-SA-5-PEG(10 kDa) (7.4 mg). ST3Gal3 (1.58U) was added to both tubes and the reaction mixtures were incubated at 32°C for 96 hrs. The reaction was monitored by SDS-PAGE gel using reagents and conditions described by Invitrogen. When the reaction was complete, the reaction mixture was purified using a Toso Haas TSK-Gel-3000 preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The combined fractions containing the product were concentrated at 4°C in Centricon-Plus-20 centrifugal filters (Millipore, Bedford, MA) and the concentrated solution reformulated to yield 1.97 mg (bicinchoninic acid protein assay, BCA assay, Sigma-Aldrich, St. Louis MO) of Factor VIIa-PEG. The product of the reaction was analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples were dialyzed against water and analyzed by MALDI-TOF. Figure 158 shows the MALDI results for native Factor VIIa. Figure 159 contains the MALDI results for Factor VIIa PEGylated with 1 kDa PEG where peak of Factor VIIa PEGylated with 1KDa PEG is evident. Figure 160 contains the MALDI results for Factor VIIa PEGylated with 10 kDa PEG where a peak for Factor VIIa PEGylated with 10 kDa PEG is evident. Figure 161 depicts the SDS-PAGE analysis of all of the reaction products, where a band for Factor VIIa-SA-PEG (10 kDa) is evident.

Follicle Stimulating Hormone (FSH)

23. GlycoPEGylation of human pituitary-derived FSH

This example illustrates the assembly of a conjugate of the invention. Follicle Stimulating Hormone (FSH) is desialylated and then conjugated with CMP-(sialic acid)-PEG.

Desialylation of Follicle Stimulating Hormone. Follicle Stimulating Hormone (FSH) (Human Pituitary, Calbiochem Cat No. 869001), 1 mg, was dissolved in 500 µL 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, 5 mM CaCl₂. This solution, 375 µL, was transferred to a small plastic tube and to it was added 263 mU Neuraminidase II (*Vibrio cholerae*). The reaction mixture was shaken gently for 15 hours at 32 °C. The reaction mixture was added to

N-(*p*-aminophenyl)oxamic acid-agarose conjugate, 600 μ L, pre-equilibrated with 50 mM Tris-HCl pH 7.4, 150 mM NaCl and 0.05% NaN₃ and gently rotated 6.5 hours at 4 °C. The suspension was centrifuged for 2 minutes at 14,000 rpm and the supernatant was collected. The beads were washed 5 times with 0.5 mL of the buffer and all supernatants were pooled.

The enzyme solution was dialyzed (7000 MWCO) for 15 hours at 4 °C with 2 L of a solution containing 50 mM Tris-HCl pH 7.4, 1 M NaCl, 0.05% NaN₃, and then twice for 4 hours at 4 °C into 50 mM Tris-HCl pH 7.4, 1 M NaCl, 0.05% NaN₃. The solution was concentrated to 2 μ g/ μ L by Speed Vac and stored at -20 °C. Reaction samples were analyzed by IEF gels (pH 3-7) (Invitrogen) (Figure 162).

Preparation of human pituitary-derived SA-FSH and PEG-SA-Follicle

Stimulating Hormone. Desialylated FSH (100 μ g, 50 μ L) and CMP-sialic acid or CMP-SA-PEG (1 kDa or 10 kDa) (0.05 μ mol) were dissolved in 13.5 μ L H₂O (adjusted to pH 8 with NaOH) in 0.5 mL plastic tubes. The tubes were vortexed briefly and 40 mU ST3Gal3 (36.5 μ L) was added (total volume 100 μ L). The tubes were vortexed again and shaken gently for 24 hours at 32 °C. The reactions were stopped by freezing at -80 °C. Reaction samples of 15 μ g were analyzed by SDS-PAGE (Figure 163), IEF gels (Figure 164) and MALDI-TOF. Native FSH was also analyzed by SDS-PAGE (Figure 165)

Analysis of SDS PAGE and IEF Gels of Reaction Products. Novex Tris-Glycine 8-16% 1 mm gels for SDS PAGE analysis were purchased from Invitrogen. 7.5 μ L (15 μ g) of FSH reaction samples were diluted with 5 μ L of 50 mM Tris-HCl, pH 7.4, 150 mM NaCl, 0.05% NaN₃ buffer, mixed with 15 μ L sample loading buffer and 1 μ L 9 M μ -mercaptoethanol and heated for 6 minutes at 85 °C. Gels were run as directed by Invitrogen and stained with Colloidal Blue Stain (Invitrogen).

FSH samples (15 μ g) were diluted with 5 μ L Tris buffer and mixed with 15 μ L sample loading buffer (Figure 162). The samples were then applied to Isoelectric Focusing Gels (pH 3-7) (Invitrogen) (Figure 165). Gels were run and fixed as directed by Invitrogen and then stained with Colloidal Blue Stain.

24. GlycoPEGylation of recombinant FSH produced recombinantly in CHO cells

This example illustrates the assembly of a conjugate of the invention. Desialylated FSH was conjugated with CMP-(sialic acid)-PEG.

- 5 **Preparation of recombinant Asialo-Follicle Stimulation Hormone.** Recombinant Follicle Stimulation Hormone (rFSH) produced from CHO was used in these studies. The 7,500 IU of rFSH was dissolved in 8 mL of water. The FSH solution was dialyzed in 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, 5 mM CaCl₂ and concentrated to 500 μ L in a Centricon Plus 20 centrifugal filter. A portion of this solution (400 μ L) (~0.8 mg FSH) was transferred to a
- 10 small plastic tube and to it was added 275 mU Neuraminidase II (*Vibrio cholerae*). The reaction mixture was mixed for 16 hours at 32 °C. The reaction mixture was added to prewashed N-(p-aminophenyl)oxamic acid-agarose conjugate (800 μ L) and gently rotated for 24 hours at 4 °C. The mixture was centrifuged at 10,000 rpm and the supernatant was collected. The beads were washed 3 times with 0.6 mL Tris-EDTA buffer, once with 0.4 mL
- 15 Tris-EDTA buffer and once with 0.2 mL of the Tris-EDTA buffer and all supernatants were pooled. The supernatant was dialyzed at 4 °C against 2 L of 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN₃ and then twice more against 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN₃. The dialyzed solution was then concentrated to 420 μ L in a Centricon Plus 20 centrifugal filter and stored at -20 °C.
- 20 Native and desialylated rFSH samples were analyzed by SDS-PAGE and IEF (Figure 166). Novex Tris-Glycine 8-16% 1 mm gels were purchased from Invitrogen. Samples (7.5 μ L, 15 μ g) samples were diluted with 5 μ L of 50 mM Tris-HCl, pH 7.4, 150 mM NaCl, 0.05% NaN₃ buffer, mixed with 15 μ L sample loading buffer and 1 μ L 9 M β -mercaptoethanol and heated for 6 minutes at 85 °C. Gels were run as directed by Invitrogen
- 25 and stained with Colloidal Blue Stain (Invitrogen). Isoelectric Focusing Gels (pH 3-7) were purchased from Invitrogen. Samples (7.5 μ L, 15 μ g) were diluted with 5 μ L Tris buffer and mixed with 15 μ L sample loading buffer. Gels were loaded, run and fixed as directed by Invitrogen. Gels were stained with Colloidal Blue Stain. Samples of native and desialylated FSH were also dialyzed against water and analyzed by MALDI-TOF.
- 30 **Sialyl-PEGylation of recombinant Follicle Stimulation Hormone.** Desialylated FSH (100 μ g, 54 μ L) and CMP-SA-PEG (1 kDa or 10 kDa) (0.05 μ mol) were dissolved in 28

μL 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN₃, pH 7.2 in 0.5 mL plastic tubes. The tubes were vortexed briefly and 20 mU of ST3Gal3 was added (total volume 100 μL). The tubes were vortexed again, mixed gently for 24 hours at 32 °C and the reactions stopped by freezing at -80 °C. Samples of this reaction were analyzed as described above by SDS-PAGE gels

5 (Figure 167), IEF gels (Figure 168) and MALDI-TOF MS.

MALDI was also performed on the PEGylated rFSH. During ionization, SA-PEG is eliminated from the N-glycan structure of the glycoprotein. Native FSH gave a peak at 13928; AS-rFSH (13282); resialylated r-FSH (13332); PEG1000-rFSH (13515; 14960 (1); 16455 (2); 17796 (3); 19321 (4)); and PEG 10000 (23560 (1); 34790 (2); 45670 (3); and 10 56760 (4)).

25. Pharmacokinetic Study of GlycoPEGylated FSH

This example sets forth the *in vivo* testing of the pharmacokinetic properties glycoPEGylated Follicle Stimulating Hormone (FSH) prepared according to the methods of 15 the invention as compared to non-PEGylated FSH.

FSH, FSH-SA-PEG (1 kDa) and FSH-SA-PEG (10 kDa) were radioiodinated using standard conditions (Amersham Biosciences, Arlington Heights, IL) and formulated in phosphate buffered saline containing 0.1% BSA. After dilution in phosphate buffer to the appropriate concentration, each of the test FSH proteins (0.4 μg, each) was injected 20 intravenously into female Sprague Dawley rats (250-300 g body weight) and blood drawn at time points from 0 to 80 hours. Radioactivity in blood samples was analyzed using a gamma counter and the pharmacokinetics analyzed using standard methods (Figure 169). FSH was cleared from the blood much more quickly than FSH-PEG(1 kDa), which in turn was clear somewhat more quickly than FSH-PEG(10 kDa).

25

26. Sertoli Cell Bioassay for *In Vitro* Activity of GlycoPEGylated FSH

This example sets forth a bioassay for follicle stimulating hormone (FSH) activity based on cultured Sertoli cells. This assay is useful to determine the bioactivity of FSH after glycan remodeling, including glycoconjugation.

30 This bioassay is based on the dose-response relationship that exists between the amount of estradiol produced when FSH, but not lutenizing hormone (LH), is added to

cultured Sertoli cells obtained from immature old rats. Exogenous testosterone is converted to 17 β -estradiol in the presence of FSH.

Seven to 10 days old Sprague-Dawley rats were used to obtain Sertoli cells. After sacrifice, testes were decapsulated and tissue was dispersed by incubation in collagenase (1 mg/ml), trypsin (1mg/ml), hyaluronidase (1 mg/ml) and DNases (5 μ g/ml) for 5 to 10 min. The tubule fragments settled to the bottom of the flask and were washed in PBS (1x). The tubule fragments were reincubated for 20 min with a media containing the same enzymes: collagenase (1 mg/ml), trypsin (1mg/ml), hyaluronidase (1 mg/ml) and DNases (5 μ g/ml).

The tubule fragments were homogenized and plated into a 24 well plate in a serum free media. 5 x 10⁵ cells were dispersed per well. After 48h incubation at 37° C and 5% CO₂, fresh media was added to the cells. Composition of the serum free media: DMEM (1 vol), Ham's F10 nutrient mixture (1 vol), insulin 1 μ g/ml, Transferrin 5 μ g/ml, EGF 10 ng/ml, T4 20 pg/ml, Hydrocortisone 10⁻⁸ M, Retinoic acid 10⁻⁶ M.

The stimulation experiment consists of a 24 hour incubation with standard FSH or samples at 37°C and 5% CO₂. The mean intra-assay coefficient of variation is 9% and the mean inter-assay coefficient of variation is 11%.

The 17B-estradiol Elisa Kit DE2000 (R&D Systems, Minneapolis, MN) was used to quantify the level of estradiol after incubation with FSH, FSH-SA-PEG (1 kDa) and FSH-SA-PEG (10 kDa).

The procedure was as follows: 100 μ l of Estradiol Standard (provided with kit and prepared as per instructions with kit) or sample was pipetted into wells of 17B-estradiol Elisa plate(s); 50 μ l of 17B-estradiol Conjugate (provided with kit, prepared as per instructions with kit) was added to each well; 50 μ l of 17B-estradiol antibody solution (provided with kit and prepared as per instructions with kit) was added to each well; plates were incubated for 2 hour at room temperature at 200 rpm; the liquid was aspirated from each well; the wells were washed 4 times using the washing solution; all the liquid was removed from the wells; 200 μ l of pNPP Substrate (provided with kit and prepared as per instructions with kit) was added to all wells and incubated for 45 min; 50 μ l of Stop solution (provided with kit and prepared as per instructions with kit) was added and the plates were read it at 405 nm (Figure 170).

While FSH-PEG(10 kDa) exhibited a modest stimulation of Sertoli cells, at 1 μ g/ml, FSH-PEG(1 kDa) stimulated Sertoli cells up to 50% more than unPEGylated FSH.

27. Steelman-Pohley Bioassay of *In Vivo* Activity of GlycoPEGylated FSH

In this example, the Steelman-Pohley bioassay (Steelman and Pohley, 1953, Endocrinology 53:604-615) was used to determine the *in vivo* activity of glycoPEGylated FSH. The Steelman-Pohley assay uses the change in ovary weight of a rat to measure the *in vivo* activity of FSH that is coinjected with human chorionic gonadotropin.

The Steelman-Pohley bioassay was performed according to the protocol described in Christin-Maitre et al. (2000, Methods 21:51-57). Seventy female Sprague-Dawley Rats (Charles River Laboratories, Wilmington, MA), aged 21 to 22 days, were housed in the testing facility for at least 5 days before the beginning of the assay procedure. Throughout the procedure, the animal room was climate controlled at 18 to 26°C, 30 to 70% relative humidity, and 12 hr. artificial light/12 hr. dark. All animals were fed Certified Rodent Chow (Harlan Teklad, Madison WI) or the equivalent, and water, both *ad libitum*. Animal procedures were performed at Calvert Preclinical Services, Inc. (Olyphant, PA).

Recombinant FSH was expressed in CHO cells, purified by standard techniques and glycoPEGylated with PEG (1 kDa). The rats were divided into seven test groups, with ten animals per group. On days -1 and 0, animals of all groups were subcutaneously injected with 20 I.U. of human chorionic gonadotropin (HCG) in 0.5 ml of 0.9 % NaCl. On days 1, 2 and 3, the control animals were subcutaneously injected with a dose of 0.5 ml containing 20 I.U. HCG in 0.9% NaCl, while in the other groups, the HCG dose was augmented with either rFSH or rFSH-SA-PEG (1 kDa) at either 0.14 µg, 0.4 µg or 1.2 µg per dose. On day 4, the animals were euthanized by CO₂ inhalation. The ovaries were removed, trimmed and weighed. The average ovary weight was determined for each group.

Figure 171 presents the average ovary weight of the test groups on day 4. The groups receiving HCG alone (control) or the low dose (0.14 µg) of either rFSH or rFSH-SA-PEG (1 kDa) had ovary weights that were roughly equivalent. The groups receiving the medium (0.4 µg) or high (1.2 µg) doses of rFSH or rFSH-SA-PEG (1 kDa) had ovary weights roughly twice that of the control group. At the medium dose (0.4 µg), the glycoPEGylated rFSH had roughly the same *in vivo* activity (as determined by ovary weight) as the unPEGylated rFSH.

At the high dose (1.2 µg), the glycoPEGylated rFSH had somewhat higher *in vivo* activity than the unPEGylated rFSH.

G-CSF

28. GlycoPEGylation of G-CSF produced in CHO cells

Preparation of Asialo-Granulocyte-Colony Stimulation Factor (G-CSF). G-CSF produced in CHO cells is dissolved at 2.5 mg/mL in 50 mM Tris 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, 5 mM CaCl₂ and concentrated to 500 µL in a Centricon Plus 20 centrifugal filter. The solution is incubated with 300 mU/mL Neuraminidase II (*Vibrio cholerae*) for 16 hours at 32 °C. To monitor the reaction a small aliquot of the reaction is diluted with the appropriate buffer and a IEF gel performed. The reaction mixture is then added to prewashed N-(*p*-aminophenyl)oxamic acid-agarose conjugate (800 µL/mL reaction volume) and the washed beads gently rotated for 24 hours at 4 °C. The mixture is centrifuged at 10,000 rpm and the supernatant was collected. The beads are washed 3 times with Tris-EDTA buffer, once with 0.4 mL Tris-EDTA buffer and once with 0.2 mL of the Tris-EDTA buffer and all supernatants are pooled. The supernatant is dialyzed at 4 °C against 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN₃ and then twice more against 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN₃. The dialyzed solution is then concentrated using a Centricon Plus 20 centrifugal filter and stored at -20 °C. The conditions for the IEF gel were run according to the procedures and reagents provided by Invitrogen. Samples of native and desialylated G-CSF are dialyzed against water and analyzed by MALDI-TOF MS.

Preparation of G-CSF-(alpha2,3)-Sialyl-PEG. Desialylated G-CSF was dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN₃, pH 7.2. The solution is incubated with 1 mM CMP-sialic acid-PEG and 0.1 U/mL of ST3Gal1 at 32°C for 2 days. To monitor the incorporation of sialic acid-PEG, a small aliquot of the reaction had CMP-SA-PEG-fluorescent ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the peptide is quantitated using an in-line fluorescent detector. After 2 days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis

according to the procedures and reagents supplied by Invitrogen. Samples of native and PEGylated G-CSF are dialyzed against water and analyzed by MALDI-TOF MS.

Preparation of G-CSF-(alpha2,8)-Sialyl-PEG. G-CSF produced in CHO cells, which contains an alpha2,3-sialylated O-linked glycan, is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN₃, pH 7.2. The solution is incubated with 1 mM CMP-sialic acid-PEG and 0.1 U/mL of CST-II at 32°C for 2 days. To monitor the incorporation of sialic acid-PEG, a small aliquot of the reaction has CMP-SA-PEG-fluorescent ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the peptide is quantitated using an in-line fluorescent detector. After 2 days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples of native and PEGylated G-CSF are dialyzed against water and analyzed by MALDI-TOF MS.

Preparation of G-CSF-(alpha2,6)-Sialyl-PEG. G-CSF, containing only O-linked GalNAc, is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN₃, pH 7.2. The solution is incubated with 1 mM CMP-sialic acid-PEG and 0.1 U/mL of ST6GalNAcI or II at 32°C for 2 days. To monitor the incorporation of sialic acid-PEG, a small aliquot of the reaction has CMP-SA-PEG-fluorescent ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the peptide is quantitated using an in-line fluorescent detector. After 2 days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples of native and PEGylated G-CSF are dialyzed against water and analyzed by MALDI-TOF MS.

G-CSF produced in CHO cells was treated with Arthrobacter sialidase and was then purified by size exclusion on Superdex75 and was treated with ST3GalI or ST3Gal2 and then with CMP-SA-PEG 20Kda. The resulting molecule was purified by ion exchange and

gel filtration and analysis by SDS PAGE demonstrated that the PEGylation was complete. This is the first demonstration of glycoPEGylation of an O-linked glycan.

Glucocerebrosidase

5 29. Glucocerebrosidase-mannose-6-phosphate produced in CHO cells

This example sets forth the procedure to glycoconjugate mannose-6-phosphate to a peptide produced in CHO cells such as glucocerebrosidase.

- Preparation of asialo-glucoceramide.** Glucocerebrosidase produced in CHO cells is dissolved at 2.5 mg/mL in 50 mM Tris 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, and is
- 10 incubated with 300 mU/mL sialidase-agarose conjugate for 16 hours at 32 °C. To monitor the reaction a small aliquot of the reaction is diluted with the appropriate buffer and a IEF gel and SDS-PAGE performed according to Invitrogen procedures. The mixture is centrifuged at 10,000 rpm and the supernatant is collected. The beads are washed 3 times with Tris-EDTA buffer, once with 0.4 mL Tris-EDTA buffer, and once with 0.2 mL of the Tris-EDTA buffer.
- 15 All supernatants are pooled. The supernatant is dialyzed at 4 °C against 50 mM Tris-HCl pH 7.4, 1 M NaCl, 0.05% NaN₃ and then twice more against 50 mM Tris-HCl pH 7.4, 1 M NaCl, 0.05% NaN₃. The dialyzed solution is then concentrated using a Centricon Plus 20 centrifugal filter. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed
- 20 against water and analyzed by MALDI-TOF MS.

Preparation of Glucocerebrosidase-SA-linker-Mannose-6-phosphate (procedure

- 1).** Asialo-glucocerebrosidase from above is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN₃, pH 7.2. The solution is incubated with 1 mM CMP-sialic acid-linker-Man-6-phosphate and 0.1 U/mL of ST3Gal3 at 32°C for 2 days. To monitor the
- 25 incorporation of sialic acid-linker-Man-6-phosphate, a small aliquot of the reaction had CMP-SA-PEG-fluorescent ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas TSK-Gel-3000 analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the peptide is quantitated using an in-line fluorescent detector. When the reaction is complete, the reaction mixture is purified
- 30 using a Toso Haas TSK-Gel-3000 preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using

SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

Preparation of Glucocerebrosidase-SA-linker-Mannose-6-phosphate (procedure

- 2). Glucocerebrosidase, produced in CHO but incompletely sialylated, is dissolved at 2.5
5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN₃, pH 7.2. The solution is incubated
with 1 mM CMP-sialic acid-linker-Man-6-phosphate and 0.1 U/mL of ST3Gal3 at 32°C for 2
days. To monitor the incorporation of sialic acid-linker-Man-6-phosphate, a small aliquot of
the reaction had CMP-SA-PEG-fluorescent ligand added; the label incorporated into the
peptide is separated from the free label by gel filtration on a Toso Haas TSK-Gel-3000
10 analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the
peptide is quantitated using an in-line fluorescent detector. When the reaction is complete,
the reaction mixture is purified using a Toso Haas TSK-Gel-3000 preparative column using
PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the
reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and
15 reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by
MALDI-TOF MS.

30. Glucocerebrosidase-transferrin

- This example sets forth the procedures for the glycoconjugation of proteins, and in
20 particular, transferrin is glycoconjugated to glucocerebrosidase. The GlcNAc-ASN structures
are created on glucoceraminidase, and Transferrin-SA-Linker-Gal-UDP is conjugated to
GNDF GlcNAc-ASN structures using galactosyltransferase.

- Preparation of GlcNAc-glucocerebrosidase (Cerezyme™).** Cerezyme™
(glucocerebrosidase) produced in CHO cells is dissolved at 2.5 mg/mL in 50 mM Tris 50 mM
25 Tris-HCl pH 7.4, 0.15 M NaCl, and is incubated with 300 mU/mL Endo-H-agarose conjugate
for 16 hours at 32 °C. To monitor the reaction a small aliquot of the reaction is diluted with
the appropriate buffer and a IEF gel and SDS-PAGE performed according to Invitrogen
procedures. The mixture is centrifuged at 10,000 rpm and the supernatant is collected. The
beads are washed 3 times with Tris-EDTA buffer, once with 0.4 mL Tris-EDTA buffer and
30 once with 0.2 mL of the Tris-EDTA buffer and all supernatants are pooled. The supernatant
is dialyzed at 4 °C against 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN₃ and then twice

more against 50 mM Tris-HCl pH 7.4, 1 M NaCl, 0.05% NaN₃. The dialyzed solution is then concentrated using a Centricon Plus 20 centrifugal filter. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF

5 MS.

Preparation of Transferrin-SA-Linker-Gal-glucocerebrosidase. Transferrin-SA-Linker-Gal-UDP from above is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 5 mM MnCl₂, 0.05% NaN₃, pH 7.2. The solution is incubated with 2.5 mg/mL GlcNAc-glucocerebrosidase and 0.1 U/mL of galactosyltransferase at 32°C for 2 days. To monitor the incorporation of glucocerebrosidase, the peptide is separated by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1) and the product detected by UV absorption. The reaction mixture is then purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

GM-CSF

31. Generation and PEGylation of GlcNAc-ASN Structures: GM-CSF produced in *Saccharomyces*

20

This example sets forth the preparation of Tissue-type Activator with PEGylated GlcNAc-Asn structures.

Recombinant GM-CSF expressed in yeast is expected to contain 2 N-linked and 2 O-linked glycans. The N-linked glycans should be of the branched mannan type. This recombinant glycoprotein is treated with an endoglycosidase from the group consisting of endoglycosidase H, endoglycosidase-F1, endoglycosidase-F2, endoglycosidase-F3, endoglycosidase-M either alone or in combination with mannosidases I, II and III to generate GlcNAc nubs on the asparagine (Asn) residues on the peptide/protein backbone.

The GlcNAc-Asn structures on the peptide/protein backbone is then be modified with galactose or galactose-PEG using UDP-galactose or UDP-galactose-6-PEG, respectively, and a galactosyltransferase such as GalT1. In one case the galactose-PEG is the terminal residue.

30

In the second case the galactose is further modified with SA-PEG using a CMP-SA-PEG donor and a sialyltransferase such as ST3GalIII. In another embodiment the GlcNAc-Asn structures on the peptide/protein backbone can be galactosylated and sialylated as described above, and then further sialylated using CMP-SA-PEG and an α 2,8-sialyltransferase such as the enzyme encoded by the *Campylobacter jejuni* cst-II gene.

Herceptin™

32. Glycoconjugation of mithramycin to Herceptin™

This example sets forth the procedures to glycoconjugate a small molecule, such as mithramycin to Fc region glycans of an antibody molecule produced in mammalian cells. Here, the antibody Herceptin™ is used, but one of skill in the art will appreciate that the method can be used with many other antibodies.

Preparation of Herceptin™-Gal-linker-mithramycin. Herceptin™ is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 5 mM MnCl₂, 0.05% NaN₃, pH 7.2. The solution is incubated with 1 mM UDP-galactose-linker-mithramycin and 0.1 U/mL of galactosyltransferase at 32°C for 2 days to introduce the mithramycin in the Fc region glycans. To monitor the incorporation of galactose, a small aliquot of the reaction has ¹⁴C-galactose-UDP ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The radioactive label incorporation into the peptide is quantitated using an in-line radiation detector.

When the reaction is complete, the reaction mixture is purified using a Toso Haas TSK-Gel-3000 preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The fractions containing product are combined, concentrated, buffer exchanged and then freeze-dried. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

Interferon α and Interferon β 33. GlycoPEGylation of Proteins expressed in Mammalian or Insect Systems:
EPO, Interferon α and Interferon β

5 This example sets forth the preparation of PEGylated peptides that are expressed in mammalian and insect systems.

Preparation of acceptor from mammalian expression systems. The peptides to be glycoPEGylated using CMP-sialic acid PEG need to have glycans terminating in galactose. Most peptides from mammalian expression systems will have terminal sialic acid that first
10 needs to be removed.

Sialidase digestion. The peptide is desialylated using a sialidase. A typical procedure involves incubating a 1 mg/mL solution of the peptide in Tris-buffered saline, pH 7.2, with 5 mM CaCl_2 added, with 0.2 U/mL immobilized sialidase from *Vibrio cholera* (Calbiochem) at 32°C for 24 hours. Microbial growth can be halted either by sterile filtration
15 or the inclusion of 0.02% sodium azide. The resin is then removed by centrifugation or filtration, and then washed to recover entrapped peptide. At this point, EDTA may be added to the solution to inhibit any sialidase that has leached from the resin.

Preparation from insect expression systems. EPO, interferon-alpha, and interferon-beta may also be expressed in non-mammalian systems such as yeast, plants, or
20 insect cells. The peptides to be glycoPEGylated using CMP-sialic acid PEG need to have glycans terminating in galactose. The majority of the N-glycans on peptides expressed in insect cells, for example, are the trimannosyl core. These glycans are first built out to glycans terminating in galactose before they are acceptors for sialyltransferase.

Building acceptor glycans from trimannosyl core. Peptide (1 mg/mL) in Tris-buffered saline, pH 7.2, containing 5 mM MnCl_2 , 5 mM UDP-glcNAc, 0.05 U/mL GLCNACT I, 0.05 U/mL GLCNACT II, is incubated at 32°C for 24 hours or until the reaction is substantially complete. Microbial growth can be halted either by sterile filtration
25 or the inclusion of 0.02% sodium azide. After buffer exchange to remove UDP and other small molecules, UDP-galactose and MnCl_2 are each added to 5 mM, galactosyltransferase is added to 0.05 U/mL, and is incubated at 32°C for 24H or until the reaction is substantially
30

complete. Microbial growth can be halted either by sterile filtration or the inclusion of 0.02% sodium azide. The peptides are then ready for glycoPEGylation.

Building O-linked glycans. A similar strategy may be employed for interferon alpha to produce enzymatically the desired O-glycan Gal-GalNAc. If necessary, GalNAc linked to serine or threonine can be added to the peptide using appropriate peptide GalNAc transferases (e.g. GalNAc T1, GalNAc T2, T3, T4, etc.) and UDP-GalNAc. Also, if needed, galactose can be added using galactosyltransferase and UDP-galactose.

GlycoPEGylation using sialyltransferase. The glycopeptides (1 mg/mL) bearing terminal galactose in Tris buffered saline + 0.02% sodium azide are incubated with CMP-SA-PEG (0.75 mM) and 0.4 U/mL sialyltransferase (ST3Gal3 or ST3Gal4 for N-glycans on EPO and interferon beta; ST3Gal4, or ST3Gal1 for O-glycans on interferon alpha) at 32°C for 24 hours. Other transferases that may work include the 2,6 sialyltransferase from *Photobacterium damsella*. The acceptor peptide concentration is most preferably in the range of 0.1 mg/mL up to the solubility limit of the peptide. The concentration of CMP-SA-PEG should be sufficient for there to be excess over the available sites, but not so high as to cause peptide solubility problems due to the PEG, and may range from 50 µM up to 5 mM, and the temperature may range from 2°C up to 40°C. The time required for complete reaction will depend on the temperature, the relative amounts of enzyme to acceptor substrate, the donor substrate concentration, and the pH.

34. GlycoPEGylation of Interferon α produced in CHO cells

Preparation of Asialo-Interferon α . Interferon alpha produced from CHO cells is dissolved at 2.5 mg/mL in 50 mM Tris 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, 5 mM CaCl₂ and concentrated to 500 µL in a Centricon Plus 20 centrifugal filter. The solution is incubated with 300 mU/mL Neuraminidase II (*Vibrio cholerae*) for 16 hours at 32 °C. To monitor the reaction a small aliquot of the reaction is diluted with the appropriate buffer and a IEF gel performed. The reaction mixture is then added to prewashed N-(*p*-aminophenyl)oxamic acid-agarose conjugate (800 µL/mL reaction volume) and the washed beads gently rotated for 24 hours at 4 °C. The mixture is centrifuged at 10,000 rpm and the supernatant was collected. The beads are washed 3 times with Tris-EDTA buffer, once with 0.4 mL Tris-EDTA buffer and once with 0.2 mL of the Tris-EDTA buffer and all

supernatants were pooled. The supernatant is dialyzed at 4 °C against 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN₃ and then twice more against 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN₃. The dialyzed solution is then concentrated using a Centricon Plus 20 centrifugal filter and stored at -20 °C. The conditions for the IEF gel are run according to the procedures and reagents provided by Invitrogen. Samples of native and desialylated G-CSF are dialyzed against water and analyzed by MALDI-TOF MS.

Preparation of Interferon-alpha-(alpha2,3)-Sialyl-PEG. Desialylated interferon-alpha is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN₃, pH 7.2. The solution is incubated with 1 mM CMP-sialic acid-PEG and 0.1 U/mL of ST3Gal1 at 32°C for 2 days. To monitor the incorporation of sialic acid-PEG, a small aliquot of the reaction had CMP-SA-PEG-fluorescent ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the peptide is quantitated using an in-line fluorescent detector. After 2 days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples of native and desialylated Interferon-alpha are dialyzed against water and analyzed by MALDI-TOF MS.

Preparation of Interferon-alpha-(alpha2,8)-Sialyl-PEG. Interferon-alpha produced in CHO, which contains an alpha2,3-sialylated O-linked glycan, is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN₃, pH 7.2. The solution is incubated with 1 mM CMP-sialic acid-PEG and 0.1 U/mL of CST-II at 32°C for 2 days. To monitor the incorporation of sialic acid-PEG, a small aliquot of the reaction has CMP-SA-PEG-fluorescent ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the peptide is quantitated using an in-line fluorescent detector. After 2 days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis

according to the procedures and reagents supplied by Invitrogen. Samples of native and PEGylated interferon-alpha are dialyzed against water and analyzed by MALDI-TOF MS.

Preparation of Interferon-alpha-(alpha2,6)-Sialyl-PEG. Interferon-alpha, containing only O-linked GalNAc, was dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN₃, pH 7.2. The solution is incubated with 1 mM CMP-sialic acid-PEG and 0.1 U/mL of ST6GalNAcI or II at 32°C for 2 days. To monitor the incorporation of sialic acid-PEG, a small aliquot of the reaction had CMP-SA-PEG-fluorescent ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the peptide is quantitated using an in-line fluorescent detector. After 2 days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples of native and PEGylated interferon-alpha are dialyzed against water and analyzed by MALDI-TOF MS.

35. GlycoPEGylation of Interferon-β-1a with PEG (10 kDa) and PEG (20 kDa)

This example illustrates a procedure PEGylate Interferon-β with either PEG (10 kDa) or PEG (20 kDa).

Briefly, Interferon-β-1a (INF-β) was obtained from Biogen (Avonex™). The INF-β was first purified by Superdex-75 chromatography. The INF-β was then desialylated with *Vibrio cholerae* sialidase. The INF-β was then PEGylated with SA-PEG (10 kDa) or SA-PEG (20 kDa) and purified with Superdex-200 chromatography.

Superdex-75 chromatography purification. INF-β (150 μg) was applied to a Superdex-75 column (Amersham Biosciences, Arlington Heights, IL) and eluted with PBS with 0.5 M NaCl, 0.02 Tween-20, 20 mM histidine and 10% glycerol. The eluant was monitored for absorbance at 280 nm (Figure 172A and 172B) and fractions were collected. Peaks 4 and 5 were pooled, concentrated in an Amicon Ultra 15 spin filter (Millipore, Billerica, MA), and the buffer was exchanged to TBS with 5 mM CaCl₂, 0.02% Tween-20, 20 mM histidine and 10% glycerol.

Sialidase Reaction. The INF- β was then desialylated with *Vibrio cholera* sialidase (70 mU/ml, CALBIOCHEM®, EMD Biosciences, Inc., San Diego, CA) on agarose in TBS with 5 mM CaCl₂, 0.02% Tween-20, 20 mM histidine and 10% glycerol. The reaction was carried out at 32°C for 18 hours. The INF- β was removed from the agarose with a 0.22 μ m Spin-X™ filter (Corning Technology, Inc., Norcross, GA). Figure 173A depicts the MALDI analysis of glycans released from native INF- β . The native INF- β has many glycoforms containing terminal sialic acid moieties. Figure 173B depicts the MALDI analysis of glycans released from desialylated INF- β . The desialylated INF- β has primarily one glycoform which is bi-antennary with terminal galactose moieties.

Lectin Dot-Blot Analysis of Sialylation. Samples of the INF- β from the desialidase reaction were dot-blotted onto nitrocellulose and then blocked with Tris buffered saline (TBS: 0.05M Tris, 0.15M NaCl, pH 7.5) and DIG kit (glycan differentiation kit available from Roche #1 210 238) blocking buffer. Some of the blots were incubated with *Maackia amurensis* agglutinin (MAA) labeled with digoxigenin (DIG) (Roche Applied Science, Indianapolis, IL) to detect α 2,3-sialylation of INF- β . These blots were washed with TBS then incubated with anti-digoxigenin antibody labeled with alkaline phosphatase, then washed again with TBS and developed with NBT/X-phosphate solution, wherein NBT is 4-nitro blue tetrazolium chloride and X-phosphate is 5-bromo-4-chloro-3-indoyl phosphate. The left side of Figure 174 depicts the results of the MAA blot of INF- β after the desialylation reaction. The INF- β is partially desialylated, as indicated by the decrease in dot development as compared to native INF- β in the desialylated samples.

Other blots were incubated with *Erthrina cristagalli* lectin (ECL) labeled with biotin (Vector Laboratories, Burlingame, CA) to detect exposed galactose residues on INF- β . After incubation with 2.5 μ g/ml ECL, the blots were washed in TBS and incubated with streptavidin labeled with alkaline phosphatase. The blots were then washed again and developed. The right side of Figure 174 depicts the ECL blot after development. The increased intensity of the dot of desialylated INF- β as compared to the native INF- β indicate more exposed galactose moieties and therefore extensive desialylation.

PEGylation of Desialylated INF- β with SA-PEG (10 kDa). Desialylated INF- β (0.05 mg/ml) was PEGylated with ST3Gal3 (50 mU/ml) and CMP-SA-PEG (10 kDa) (250

μM) in an appropriate buffer of TBS + 5 mM CaCl₂, 0.02% Tween 20, 20 mM histidine, 10% glycerol for 50 hours at 32°C. Figure 175 depicts the SDS-PAGE analysis of the reaction products showing PEGylated INF-β at approximately 98 kDa.

PEGylation of Desialylated INF-β with SA-PEG (20 kDa). Desialylated INF-β (0.5 mg/ml) was PEGylated with ST3Gal3 (170 mU/ml) and CMP-SA-PEG (20 kDa) in an appropriate buffer of TBS + 5 mM CaCl₂, 0.02% Tween 20, 20 mM histidine, 10% glycerol for 50 hours at 32°C. Figure 176 depicts the SDS-PAGE analysis the products of the PEGylation reaction. The PEGylated INF-β has many higher molecular weight bands not found in the unmodified INF-β indicating extensive PEGylation.

Superdex-200 Purification of INF-β PEGylated with PEG (10 kDa). The products of the PEGylation reaction were separated on a Superdex-200 column (Amersham Biosciences, Arlington Heights, IL) in PBS with 0.5 NaCl, 0.02 Tween-20, 20 mM histidine and 10% glycerol at 1ml/min and 30 cm/hr flow. The eluant was monitored for absorbance at 280 nm (Figure 177) and fractions were collected. Peaks 3 and 4 were pooled and concentrated in an Amicon Ultra 15 spin filter.

Bioassay of INF-β PEGylated with PEG (10 kDa).

The test is inhibition of the proliferation of the lung carcinoma cell line, A549. The A549 cell line are lung carcinoma adherent cells growing in RPMI + 10% FBS at 37°C 5% CO₂. They can be obtained from ATCC # CCL-185. Wash the cells with 10 ml of PBS and remove the PBS. Add 5 ml of trypsin, incubate for 5 minutes at room temperature or 2 minutes at 37°C. When the cells are detached resuspend into 25 ml of media and count the cells. Dilute the cells at a concentration of 10000 cells/ml and add 200 ul / well (96 wells plate). Incubate for 4 hours at 37°C 5% CO₂. Prepare 1 ml of IFN B at a concentration of 0.1 ug/ml. Filter it under the hood with a 0.2 um filter. Add 100 ul per well (8 replicates = 1 lane). Incubate for 3 days (do not let the cells go to confluence). Remove 200 ul of media (only 100ul per well left). Add 25 μl of MTT (Sigma) (5 mg/ml filtered 0.22μm). Incubate for 4 hours at 37°C and 5% CO₂. Aspirate the media gently and add 100 μl of a mixture of isopropanol (100 ml and 6N HCl. Aspirate up and down to homogenize the crystal violet. Read OD 570nm (remove the background at 630 or 690 nm).

Figure 178 depicts the results of the bioassay of the peaks containing INF- β PEGylated with PEG (10 kDa) as eluted from the Superdex-200 column.

Superdex-200 Purification of INF- β PEGylated with PEG (20 kDa). The products of the PEG (20 kDa) PEGylation reaction were separated on a Superdex-200 column (Amersham Biosciences, Arlington Heights, IL) in PBS with 0.5 NaCl, 0.02 Tween-20, 20 mM histidine and 10% glycerol at 1 ml/min flow. The eluant was monitored for absorbance at 280 nm (Figure 179) and fractions were collected. Peak 3 contained most of the INF- β PEGylated with PEG (20 kDa).

Endotoxin test of INF- β PEGylated with PEG (20 kDa).

Limulus Lysate Test was performed, BioWhittaker # 50-647U

Table 24. Results of the endotoxin test of INF- β PEGylated with PEG (20 kDa).

	Concentration		
INF- β with PEG (20 kDa)	10 EU/ml	0.06 mg/ml	0.16 EU/ μ g
INF- β with PEG (20 kDa)	1 EU/ml	0.07 mg/ml	0.014 EU/ μ g
Native INF- β	40 EU/ml	0.1 mg/ml	0.4 EU/ μ g

Remicade™

36. GlycoPEGylation of Remicade™ antibody

This example sets forth the procedure to glycoPEGylate a recombinant antibody molecule by introducing PEG molecules to the Fc region glycans. Here Remicade™, a TNF-R: IgG Fc region fusion protein, is the exemplary peptide.

Preparation of Remicade™-Gal-PEG (10 kDa). Remicade™ is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 5 mM MnCl₂, 0.05% NaN₃, pH 7.2. The solution is incubated with 1 mM UDP-galactose-PEG (10 kDa) and 0.1 U/mL of galactosyltransferase at 32°C for 2 days to introduce the PEG in the Fc region glycans. To monitor the incorporation of galactose, a small aliquot of the reaction has ¹⁴C-galactose-UDP ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The radioactive label incorporation into the peptide is quantitated using an in-line radiation detector.

When the reaction is complete, the reaction mixture is purified using a Toso Haas TSK-Gel-3000 preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The fractions containing product are combined, concentrated, buffer exchanged and then freeze-dried. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

Rituxan™

37. Glycoconjugation of geldanamycin to Rituxan™

This example sets forth the glycoconjugation of a small molecule, such as geldanamycin, to the Fc region glycans of an antibody produced in CHO cells, such as Rituxan™. Here, the antibody Rituxan™ is used, but one of skill in the art will appreciate that the method can be used with many other antibodies.

Preparation of Rituxan™-Gal-linker-geldanamycin. Rituxan™ is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 5 mM MnCl₂, 0.05% NaN₃, pH 7.2. The solution is incubated with 1 mM UDP-galactose-linker-geldanamycin and 0.1 U/mL of galactosyltransferase at 32°C for 2 days to introduce the geldanamycin in the Fc region glycans. To monitor the incorporation of galactose, a small aliquot of the reaction has ¹⁴C-galactose-UDP ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The radioactive label incorporation into the peptide is quantitated using an in-line radiation detector.

When the reaction is complete, the reaction mixture is purified using a Toso Haas TSK-Gel-3000 preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The fractions containing product are combined, concentrated, buffer exchanged and then freeze-dried. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

Rnase38. Remodeling high mannose N-glycans to hybrid and complex N-glycans:
Bovine pancreatic RNase

This example sets forth the preparation of bovine pancreas RNase with hybrid or
5 complex N-glycans. The high mannose N-linked glycans of the RNase are enzymatically
digested and elaborated to create hybrid N-linked glycans. Additionally, the high mannose
N-linked glycans of the RNase are enzymatically digested and elaborated to create complex
N-linked glycans.

High mannose structures of *N*-linked oligosaccharides in glycopeptides can be
10 modified to hybrid or complex forms using the combination of α -mannosidases and
glycosyltransferases. This example summarizes the results in such efforts using a simple *N*-
Glycan as a model substrate.

Ribonuclease B (RNaseB) purified from bovine pancreas (Sigma) is a glycopeptide
consisting of 124 amino acid residues. It has a single potential *N*-glycosylation site modified
15 with high mannose structures. Due to its simplicity and low molecular weight (13.7 kDa to
15.5 kDa), ribonuclease B is a good candidate to demonstrate the feasibility of the *N*-Glycan
remodeling from high mannose structures to hybrid or complex *N*-linked oligosaccharides.
The MALDI-TOF spectrum of RNaseB (Figure 180A) and HPLC profile for the
oligosaccharides cleaved from RNaseB by *N*-Glycanase (Figure 180B) indicated that, other
20 than a small portion of the non-modified peptide, the majority of *N*-glycosylation sites of the
peptide are modified with high mannose oligosaccharides consisting of 5 to 9 mannose
residues.

Conversion of high mannose N-Glycans to hybrid N-Glycans. High mannose *N*-
Glycans were converted to hybrid *N*-Glycans using the combination of α 1,2-mannosidase,
25 GlcNAcT-I (β -1,2-*N*-acetyl glucosaminyl transferase), GalT-I (β 1,4-galactosyltransferase) and
 α 2,3-sialyltransferase /or α 2,6-sialyltransferase as shown in Figure 181.

As an example, high mannose structures in RNaseB were successfully converted to
hybrid structures.

Man₅GlcNAc₂-R was obtained from Man₅₋₉GlcNAc₂-R catalyzed by a single α 1,2-
30 mannosidase cloned from *Trichoderma reesei* (Figure 182). RNase B (1 g, about 67 μ mol)
was incubated at 30°C for 45 hr with 15 mU of the recombinant *T. reesei* α 1,2-mannosidase

in MES buffer (50 mM, pH 6.5) in a total volume of 10 mL. $\text{Man}_6\text{GlcNAc}_2$ -protein structures have been successfully converted to $\text{Man}_5\text{GlcNAc}_2$ -protein with high efficiency by the recombinant mannosidase.

Alternately, $\text{Man}_5\text{GlcNAc}_2$ -R was obtained from $\text{Man}_5\text{GlcNAc}_2$ -R catalyzed by a single α 1,2-mannosidase purified from *Aspergillus saitoi* (Figure 183). RNase B (40 μ g, about 2.7 nmol) was incubated at 37°C for 42.5 hr with 25 μ U of the commercial *A. saitoi* α 1,2-mannosidase (Glyko or CalBioChem) in NaOAc buffer (100 mM, pH 5.0) in a total volume of 20 μ L. $\text{Man}_6\text{GlcNAc}_2$ -protein structures were successfully converted to $\text{Man}_5\text{GlcNAc}_2$ -protein by the commercially available mannosidase. However, a new peak corresponding to the GlcNAc-protein appears in the spectrum, indicating the possible contamination of endoglycosidase H in the preparation. Although several mammalian α 1,2-mannosidases were required to achieve this step, the fungal α 1,2-mannosidase was very efficient to remove all α 1,2-linked mannose residues.

GlcNAcT-I then added a GlcNAc residue to the $\text{Man}_5\text{GlcNAc}_2$ -R (Figure 184). The reaction mixture after the *T. reesei* α 1,2-mannosidase reaction containing RNase B (600 μ g, about 40 nmol) was incubated with non-purified recombinant GlcNAcT-I (34 mU) in MES buffer (50 mM, pH 6.5) containing MnCl_2 (20 mM) and UDP-GlcNAc (5 mM) in a total volume of 400 μ L at 37°C for 42 hr. A GlcNAc residue was quantitatively added to $\text{Man}_5\text{GlcNAc}_2$ -protein by the recombinant GlcNAcT-I.

A Gal residue was then added using GalT 1 (Figure 185). The reaction mixture after the GnT-I reaction containing RNase B (120 μ g, about 8 nmol) was incubated at 37°C for 20 hr with 3.3 mU of the recombinant GalT-1 in Tris-HCl buffer (100 mM, pH 7.3) containing UDP-Gal (7.5 mM) and MnCl_2 (20 mM) in a total volume of 100 μ L. A Gal residue was added to about 98% of the GlcNAc- $\text{Man}_5\text{GlcNAc}_2$ -protein by the recombinant GalT 1.

The next step was the addition of a sialic acid using an α 2,3-sialyltransferase or an α 2,6-sialyltransferase (Figure 186). As an example, ST3Gal III, an α 2,3-sialyltransferase was used. The reaction mixture after the GalT-1 reaction containing RNase B (13 μ g, about 0.87 nmol) was incubated at 37°C for 16 hr with 8.9 mU of recombinant ST3Gal III in Tris-HCl buffer (100 mM, pH 7.3) containing CMP-Sialic acid (5 mM) and MnCl_2 (20 mM) in a total volume of 20 μ L. A sialic acid residue was added to about 90% of the Gal-GlcNAc-

Man₅GlcNAc₂-protein by recombinant ST3Gal III using CMP-SA as the donor. The yield can be further improved by adjusting the reaction conditions.

For convenience, no purification or dialysis step was required after each reaction described above. More interesting, GalT I and ST3Gal III can be combined in a one-pot reaction. Similar yields were obtained as compared with the separate reactions. The reaction mixture after the GlcNAcT-I reaction containing RNase B (60 µg, about 4 nmol) was incubated at 37°C for 20 hr with 1.7 mU of recombinant GalT I, 9.8 mU of recombinant ST3Gal III in Tris-HCl buffer (100 mM, pH 7.3) containing UDP-Gal (7.5 mM), CMP-sialic acid (5 mM) and MnCl₂ (20 mM) in a total volume of 60 µl.

As shown in Figure 187, SA-PEG (10 kDa) was successfully added to the RNaseB. The reaction mixture after the GalT-I reaction containing RNase B (6.7 µg, about 0.45 nmol) was dialyzed against H₂O for 1 hour at room temperature and incubated at 37°C for 15.5 hours with 55 mU of the recombinant ST3Gal III in Tris-HCl buffer (50 mM, pH 7.3) containing CMP-SA-PEG (10 kDa) (0.25 mM) and MnCl₂ (20 mM) in a total volume of 20 µl. PEG-modified sialic acid residues were successfully added to the Gal-GlcNAc-Man₅GlcNAc₂-peptide by the recombinant ST3Gal III. The yield can be further improved by adjusting the reaction conditions.

Conversion of high mannose N-Glycans to complex N-Glycans. To achieve this conversion, a GlcNAcβ1,2Man₅GlcNAc₂-peptide intermediate is obtained. As shown in Figure 188, there are at least four feasible routes to carry out the reaction from Man₅GlcNAc₂-peptide to this intermediate:

Route I: The Man₅GlcNAc₂-peptide produced by the fungal α1,2 mannosidase is a substrate of GlcNAc transferase I (GlcNAcT-I, enzyme 2) which adds one GlcNAc. The terminal α1,3- and α1,6-linked mannose residues of GlcNAcMan₅GlcNAc₂-peptide is removed by Golgi α-mannosidase II (ManII, enzyme 5). This route is a part of the natural pathway for the processing of N-linked oligosaccharides carried out in higher organisms.

Route II: Two mannose residues are first removed by an α-mannosidase (enzyme 6), then a GlcNAc is added by GlcNAcT-I (enzyme 2). Other than its natural acceptor Man₅GlcNAc₂-R, GlcNAcT-I can also recognize Man₃GlcNAc₂-R as its substrate and add one GlcNAc to the mannose core structure to form GlcNAcMan₃GlcNAc₂-peptide.

Route III: The α 1,6-linked mannose is removed by an α 1,6-mannosidase, followed by the addition of GlcNAc by GlcNAcT-I and removal of the terminal α 1,3-linked mannose by an α 1,3-mannosidase. From the experimental data obtained, GlcNAcT-I can recognize this Man₄GlcNAc₂-peptide as acceptor and add one GlcNAc residue to form

5 GlcNAcMan₄GlcNAc₂-peptide.

Route IV: Similar to Route III, α 1,3-linked mannose is removed by an α 1,3-mannosidase, followed by GlcNAcT-I reaction. Then the terminal α 1,6-linked mannose can be removed by an α 1,6-mannosidase.

After the function of GlcNAcT-I (responsible for the addition of the GlcNAc β 1,2-linked to the α 1,3-mannose on the mannose core) and GlcNAcT-II (responsible for the addition of a second GlcNAc β 1,2-linked to the α 1,6-mannose on the mannose core), the GlcNAc₂Man₃GlcNAc₂-peptide can be processed by GalT 1 and sialyltransferase to form bi-antennary complex N- Glycans. Other GlcNAc transferases such as GlcNAcT-IV, GlcNAcT-V, and/or GlcNAcT-VI (Figure 188 and Figure 189) can also glycosylate the

10 linked to the α 1,3-mannose on the mannose core) and GlcNAcT-II (responsible for the addition of a second GlcNAc β 1,2-linked to the α 1,6-mannose on the mannose core), the GlcNAc₂Man₃GlcNAc₂-peptide can be processed by GalT 1 and sialyltransferase to form bi-antennary complex N- Glycans. Other GlcNAc transferases such as GlcNAcT-IV, GlcNAcT-V, and/or GlcNAcT-VI (Figure 188 and Figure 189) can also glycosylate the

15 GlcNAc₂Man₃GlcNAc₂-peptide. Additional glycosylation by the GalT 1 and sialyltransferases will form multi-antennary complex N-glycans. The enzyme GlcNAcT-III catalyzes the insertion of a bisecting GlcNAc, thus preventing the actions of ManII and subsequent action of transferases GlcNAcT-II, GlcNAcT-IV and GlcNAcT-V.

20 Tissue-Type Plasminogen Activator (TPA)

39. Fucosylation of TPA to create Sialyl Lewis X

This example sets forth the preparation of Tissue Tissue-type Plasminogen Activator (TPA) with N-linked sialyl Lewis X antigen.

Sialylation. TPA expressed in mammalian cells will often contain a majority of the glycans terminating in sialic acid, but to ensure complete sialylation, it would be beneficial to first perform an *in vitro* sialylation. TPA in a suitable buffer (most preferably between pH 5.5 and 9, for example Tris buffered saline, pH 7.2) is incubated with CMP sialic acid and sialyltransferase for a time sufficient to convert any glycans lacking sialic acid to sialylated species. Typical conditions would be 1 mg/mL TPA, 3 mM CMP sialic acid, 0.02 U/mL

25 glycans terminating in sialic acid, but to ensure complete sialylation, it would be beneficial to first perform an *in vitro* sialylation. TPA in a suitable buffer (most preferably between pH 5.5 and 9, for example Tris buffered saline, pH 7.2) is incubated with CMP sialic acid and sialyltransferase for a time sufficient to convert any glycans lacking sialic acid to sialylated species. Typical conditions would be 1 mg/mL TPA, 3 mM CMP sialic acid, 0.02 U/mL

30 ST3Gal3, 32°C for 24 hours. Microbial growth can be halted either by sterile filtration or the

inclusion of 0.02% sodium azide. The TPA concentration is most preferably in the range 0.1 mg/mL up to the solubility limit of the peptide. The concentration of CMP-SA should be sufficient for there to be excess over the available sites, and might range from 50 μ M up to 50 mM, and the temperature from 2°C up to 40°C. The time required for complete reaction will depend on the temperature, the relative amounts of enzyme to acceptor substrate, the donor substrate concentration, and the pH. Other sialyltransferases that may be capable of adding sialic acid in 2,3 linkage include ST3Gal4; microbial transferases could also be used.

Fucosylation. Typical conditions for fucosylation would be 1 mg/mL TPA, 3 mM GDP-fucose, 0.02 U/mL FTVI, 5 mM $MnCl_2$, 32°C for 24H in Tris buffered saline.

Microbial growth can be halted either by sterile filtration or the inclusion of 0.02% sodium azide. The TPA concentration is most preferably in the range 0.1 mg/mL up to the solubility limit of the peptide. The concentration of GDP-fucose should be sufficient for there to be excess over the available sites, and might range from 50 μ M up to 50 mM, and the temperature from 2°C up to 40°C. The time required for complete reaction will depend on the temperature, the relative amounts of enzyme to acceptor substrate, the donor substrate concentration, and the pH. Other fucosyltransferases that may be capable of making sialyl Lewis x include FTVII, FTV, FTIII, as well as microbial transferases could also be used.

40. Trimming of high mannose to tri-mannose core structure: Tissue-type Plasminogen Activator produced in CHO

This example sets forth the preparation of Tissue-type Plasminogen Activator with a trimannose core by trimming back from a high mannose glycan.

Tissue-type plasminogen activator (TPA) is currently produced in Chinese Hamster Ovary (CHO) cells and contains a low amount of high mannose N-linked oligosaccharide.

The mannoses can be trimmed down using a variety of the specific mannosidases. The first step is to generate Man5GlcNAc2(Fuc0-1) from Man9GlcNAc2(Fuc0-1). This can be done using mannosidase I. Then either GlcNAcT1 (GlcNAc transferase I) is used to make GlcNAc1Man5GlcNAc2(Fuc0-1) or Mannosidase III is used to make Man3GlcNAc2(Fuc0-1). From Man3GlcNAc2(Fuc0-1), GlcNAc1Man3GlcNAc2(Fuc0-1) can be produced using GlcNAcT1 or from GlcNAc1Man5GlcNAc2(Fuc0-1), GlcNAc1Man3GlcNAc2(Fuc0-1) can be produced using Mannosidase II. GlcNAc1Man3GlcNAc2(Fuc0-1) is then converted into

GlcNAc2Man3GlcNAc2(Fuc0-1) using GlcNAcTransferase II (GlcNAcTII). The two terminal GlcNAc residues are then galactosylated using GalTI and then sialylated with SA-PEG using ST3GalIII.

- Conversely, TPA can be produce in yeast or fungal systems. Similar processing
5 would be required for fungal derived material.

41. Generation and PEGylation of GlcNAc-ASN structures: TPA produced in Yeast

- This example sets forth the preparation of PEGylated GlcNAc-Asn structures on a
10 peptide such as TPA expressed in yeast.

Yeast expression is expected to result in a TPA which contains a single N-linked mannan-type structure. This recombinant glycoprotein is first treated with endoglycosidase H to generate GlcNAc structures on the asparagine (Asn) residues on the peptide.

- The GlcNAc-Asn structures on the peptide/protein backbone are then be modified
15 with galactose or galactose-PEG using UDP-galactose or UDP-galactose-6-PEG, respectively, and a galactosyltransferase such as GalT1. In one case, the galactose-PEG is the terminal residue. In the second case, the galactose is further modified with SA-PEG using a CMP-SA-PEG donor and a sialyltransferase such as ST3GalIII. In another embodiment, the GlcNAc-Asn structures on the peptide/protein backbone may be
20 galactosylated and sialylated as described above, and then further sialylated using CMP-SA-PEG and an α 2,8-sialyltransferase such as the enzyme encoded by the *Campylobacter jejuni* cst-II gene.

Transferrin

- 25 42. GlycoPEGylation of Transferrin

This example sets forth the preparation of asialotransferrin and its sialylation with PEG-CMP-sialic acid.

- Preparation of Asialo-transferrin.** Human-derived holo-Transferrin, (10 mg) was dissolved in 500 μ L of 50 mM NaOAc, 5 mM CaCl₂, pH 5.5. To this solution was added
30 500 mU Neuraminidase II (*Vibrio cholerae*) and the reaction mixture was shaken gently for 20.5 hours at 37 °C. The reaction mixture was added to the prewashed N-(p-

aminophenyl)oxamic acid-agarose conjugate (600 μ L) and the washed beads gently rotated for 24 hours at 4 °C. The mixture was centrifuged at 10,000 rpm and the supernatant was collected. The reaction mixture was adjusted to 5 mM EDTA by addition of 100 μ L of 30 mM EDTA to the washed beads, which were gently rotated for 20 hours at 4 °C. The suspension was centrifuged for 2 minutes at 10,000 rpm and the supernatant was collected. The beads were washed 5 times with 0.35 mL of 50 mM NaOAc, 5 mM CaCl₂, 5 mM EDTA, pH 5.5 and all supernatants were pooled. The enzyme solution was dialyzed twice at 4 °C into 15 mM Tris-HCl, 1 M NaCl, pH 7.4. 0.3 mL of the transferrin solution (3.3 mL total) was removed and dialyzed twice against water. The remainder was dialyzed twice more at 4 °C against phosphate buffered saline. The dialyzed solution was stored at -20 °C. Protein samples were analyzed by IEF Electrophoresis. Samples (9 μ L, 25 μ g) were diluted with 16 μ L Tris buffer and mixed with 25 μ L of the sample loading buffer and applied to Isoelectric Focusing Gels (pH 3-7). Gels were run and fixed using standard procedures. Gels were stained with Colloidal Blue Stain.

Sialyl-PEGylation of asialo-Transferrin. Desialylated transferrin (250 μ g) and CMP-sialic acid or CMP-SA-PEG (1 kDa or 10 kDa)(0.05 μ mol) were dissolved in 69 μ L 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN₃, pH 7.2 in 1.5 mL plastic tubes. The tubes were vortexed briefly and 100 mU ST3Gal3 (90 μ L) were added (total volume 250 μ L). The tubes were vortexed again and mixed gently for 24 hours at 32 °C. The reactions were stopped by freezing at -80 °C. Novex Tris-Glycine 8-16% 1 mm gels were used for SDS PAGE analysis (Figure 190). Samples (25 μ L, 25 μ g) were mixed with 25 μ L of sample loading buffer and 0.4 μ L of β -mercaptoethanol and heated for 6 minutes at 85 °C. Gels were run using standard conditions and stained with Colloidal Blue Stain. IEF gels were also performed as described above (Figure 191). Samples were also dialyzed against water and analyzed by MALDI-TOF.

Results. MALDI was also performed. Native transferrin (78729); asialotransferrin (78197); resialylated transferrin (79626/80703); with SA-PEG 1k (79037 (1); 80961 (2); 82535 (3); 84778 (4)); with SA-PEG 5k (90003 (2); 96117 (3); 96117 (4)); with SA-PEG 10k (100336 (2); 111421 (3); 122510 (4)).

43. Transferrin-GDNF

This example sets forth the procedures for the glycoconjugation of proteins, and in particular, transferrin is glycoconjugated to GDNF. Transferrin-SA-Linker-Gal-UDP is prepared from transferrin. The galactose residue is removed from GDNF glycans, and

5 Transferrin-SA-Linker-Gal-UDP is conjugated to GDNF glycans using a galactosyltransferase.

Preparation of agalacto-GDNF. GDNF produced in NSO cells (NSO murine myeloma cells) is dissolved at 2.5 mg/mL in 50 mM Tris 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, and is incubated with 300 mU/mL beta-galactosidase-agarose conjugate for 16 hours at

10 32°C. To monitor the reaction a small aliquot of the reaction is diluted with the appropriate buffer and a IEF gel performed according to Invitrogen procedures. The mixture is centrifuged at 10,000 rpm and the supernatant is collected. The supernatant is dialyzed at 4 °C against 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN₃ and then twice more against 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN₃. The dialyzed solution is then concentrated

15 using a Centricon Plus 20 centrifugal filter and stored at -20 °C. The conditions for the IEF gel are run according to the procedures and reagents provided by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

Preparation of Transferrin-SA-Linker-Gal-UDP. Asialo-transferrin is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN₃, pH 7.2. The solution is

20 incubated with CMP-sialic acid-linker-Gal-UDP (molar amount to add 1 molar equivalent of nucleotide sugar to transferrin) and 0.1 U/mL of ST3Gal3 at 32°C for 2 days. To monitor the incorporation of sialic acid, a small aliquot of the reaction has ¹⁴C-SA-UDP ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The radioactive label

25 incorporation into the peptide is quantitated using an in-line radiation detector.

The solution is incubated with 5 mM CMP-sialic acid and 0.1 U/mL of ST3Gal3 (to cap any unreacted transferrin glycans) at 32°C for 2 days. The incorporation into the peptide is quantitated using an in-line UV detector. After 2 days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) and collecting

30 fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE

and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

- Preparation of Transferrin-SA-Linker-Gal-GDNF.** The transferrin-SA-Linker-Gal-UDP prepared as described above is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 5 mM MnCl₂, 0.05% NaN₃, pH 7.2. The solution is incubated with 2.5 mg/mL agalacto-GDNF and 0.1 U/mL of galactosyltransferase at 32°C for 2 days. To monitor the incorporation of galactose, a small aliquot of the reaction has ¹⁴C-galactose-UDP ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The radioactive label incorporation into the peptide is quantitated using an in-line radiation detector.

- When the reaction is complete, the solution is incubated with 5 mM UDP-Gal and 0.1 U/mL of galactosyltransferase (to cap any unreacted transferrin glycans) at 32°C for 2 days followed by addition of 5 mM CMP-SA and 0.1 U/mL of ST3Gal3. After 2 additional days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

- The disclosures of each and every patent, patent application, and publication cited herein are hereby incorporated herein by reference in their entirety.

- While this invention has been disclosed with reference to specific embodiments, it is apparent that other embodiments and variations of this invention may be devised by others skilled in the art without departing from the true spirit and scope of the invention. The appended claims are intended to be construed to include all such embodiments and equivalent variations.

What is claimed:

1. An EPO peptide comprising one or more glycans, having a glycoconjugate molecule covalently attached to said peptide.

2. The EPO peptide of claim 1, wherein said one or more glycans is a
5 monoantennary glycan.

3. The EPO peptide of claim 1, wherein said one or more glycans is a biantennary glycan.

4. The EPO peptide of claim 1, wherein said one or more glycans is a triantennary glycan.

10 5. The EPO peptide of claim 1, wherein said one or more glycans is at least a triantennary glycan.

6. The EPO peptide of claim 1, wherein said one or more glycans comprises at least two glycans comprising a mixture of mono or multiantennary glycans.

15 7. The EPO peptide of claim 1, wherein said one or more glycans is selected from an N-linked glycan and an O-linked glycan.

8. The EPO peptide of claim 1, wherein said one or more glycans is at least two glycans selected from an N-linked and an O-linked glycan.

9. The EPO peptide of claim 1, wherein said peptide is expressed in a cell selected from the group consisting of a prokaryotic cell and a eukaryotic cell.

20 10. The EPO peptide of claim 9, wherein said eukaryotic cell is selected from the group consisting of a mammalian cell, an insect cell and a fungal cell.

11. The EPO peptide of claim 10, wherein said fungal cell is a yeast cell.

12. A glycoPEGylated EPO peptide comprising an EPO peptide and at least one glycan and at least one poly(ethylene glycol) molecule covalently attached to said glycan,

wherein said poly(ethylene glycol) molecule is added to said EPO peptide using a glycosyltransferase.

13. The glycoPEGylated EPO peptide of claim 12, comprising at least one
5 mono-antennary glycan.

14. The glycoPEGylated EPO peptide of claim 12, wherein all of said glycans
are N-linked and are mono-antennary.

15. The glycoPEGylated EPO peptide of claim 12, wherein all of said glycans
are N-linked and at least one of said glycans comprise said poly(ethylene glycol).

16. The glycoPEGylated EPO peptide of claim 15, wherein more than one of
said glycans comprises said poly(ethylene glycol).

17. The glycoPEGylated EPO peptide of claim 12, wherein all of said glycans
are N-linked and all of said glycans comprise said poly(ethylene glycol).

18. The glycoPEGylated EPO peptide of claim 12, comprising at least three
20 mono-antennary glycans having said poly(ethylene glycol) covalently attached thereto.

19. A glycoPEGylated EPO peptide, wherein said EPO peptide comprises
three or more glycans.

20. The glycoPEGylated EPO peptide of claim 9, wherein at least one of said
glycans comprises said poly(ethylene glycol) covalently attached thereto.

21. The glycoPEGylated EPO peptide of claim 18, wherein more than one of
said glycans comprises said poly(ethylene glycol) covalently attached thereto.

22. The glycoPEGylated EPO peptide of claim 18, wherein all of said glycans comprise said poly(ethylene glycol) covalently attached thereto.

23. The glycoPEGylated EPO peptide of claim 12 wherein said poly(ethylene glycol) is linked to at least one sugar moiety selected from the group consisting of fucose (Fuc), N-acetylglucosamine (GlcNAc), galactose (Gal) and a sialic acid (SA).

24. The glycoPEGylated EPO peptide of claim 23, wherein said sialic acid is N-acetylneuraminic acid.

25. The glycoPEGylated EPO peptide of claim 12, wherein said EPO peptide does not comprise an O-linked glycan.

26. The glycoPEGylated EPO peptide of claim 12 wherein said EPO peptide comprises at least one O-linked glycan.

27. The glycoPEGylated EPO peptide of claim 26, wherein said O-linked peptide comprises said poly(ethylene glycol) covalently attached thereto.

28. The glycoPEGylated EPO peptide of claim 27, wherein said EPO peptide is recombinantly expressed in a cell.

29. The glycoPEGylated EPO peptide of claim 28, wherein said cell is selected from the group consisting of an insect cell, a fungal cell and a mammalian cell.

30. The glycoPEGylated EPO peptide of claim 29, wherein said fungal cell is a yeast cell.

31. The glycoPEGylated EPO peptide of claim 29, wherein said cell is an insect cell.

32. The glycoPEGylated EPO peptide of claim 29, wherein said cell is a yeast cell.

33. The glycoPEGylated EPO peptide of claim 29, wherein said cell is a mammalian cell.

34. The glycoPEGylated EPO peptide of claim 33, wherein said mammalian cell is a CHO cell.

35. The glycoPEGylated EPO peptide of claim 12, wherein said poly(ethylene glycol) has a molecular weight selected from the group consisting of about 1 kDa, 2 kDa, 5 kDa, 10 kDa, 20 kDa, 30 kDa and 40 kDa.

36. The glycoPEGylated EPO peptide of claim 35, wherein said poly(ethylene glycol) has a molecular weight of 20 kDa.

37. The glycoPEGylated EPO peptide of claim 12, wherein said EPO peptide is selected from the group consisting of a naturally occurring EPO peptide and a mutated EPO peptide.

38. The glycoPEGylated EPO peptide of claim 37, wherein said mutated EPO peptide comprises the amino acid sequence of SEQ ID NO:73 having at least one mutation selected from the group consisting of Arg¹³⁹ to Ala¹³⁹, Arg¹⁴³ to Ala¹⁴³ and Lys¹⁵⁴ to Ala¹⁵⁴.

39. A method of making a glycoPEGylated EPO peptide, said method comprising the step of:

(a) contacting an EPO peptide with a mixture comprising a nucleotide sugar covalently linked to poly(ethylene glycol) and a glycosyltransferase under conditions sufficient to transfer said poly(ethylene glycol) to said EPO peptide.

40. The method of claim 39, wherein the sugar of said nucleotide sugar is selected from the group consisting of fucose (Fuc), N-acetylglucosamine (GlcNAc), galactose (Gal) and a sialic acid (SA).

5 41. The method of claim 40, wherein said sialic acid is N-acetylneuraminic acid (NAN).

42. The method of claim 39, wherein said poly(ethylene glycol) has a molecular weight selected from the group consisting of about 1 kDa, 2 kDa, 5 kDa, 10 kDa,
10 20 kDa, 30 kDa and 40 kDa.

43. The method of claim 42, wherein said poly(ethylene glycol) has a molecular weight of 20 kDa.

15 44. The method of claim 39, wherein said EPO peptide is recombinantly expressed in a cell.

45. The method of claim 44, wherein said cell is selected from the group consisting of an insect cell, a fungal cell and a mammalian cell.

20 46. The method of claim 45, wherein said cell is an insect cell.

47. The method of claim 45, wherein said cell is a yeast cell.

25 48. The method of claim 45, wherein said cell is a mammalian cell.

49. The method of claim 48, wherein said mammalian cell is a CHO cell.

50. The method of claim 39, wherein said EPO peptide is selected from the
30 group consisting of a naturally occurring EPO peptide and a mutated EPO peptide.

51. The method of claim 50, wherein said mature EPO peptide has the sequence of SEQ ID NO:73.

52. The method of claim 50, wherein said mutated EPO peptide comprises the
5 amino acid sequence of SEQ ID NO: 73 having at least one mutation selected from the group consisting of Arg¹³⁹ to Ala¹³⁹, Arg¹⁴³ to Ala¹⁴³ and Lys¹⁵⁴ to Ala¹⁵⁴.

53. The method of claim 39, wherein before step (a):

(b) contacting said EPO peptide with a mixture comprising a nucleotide-N-
10 acetylglucosamine (GlcNAc) molecule and an N-acetylglucosamine transferase (GnT) for which the nucleotide-GlcNAc is a substrate under conditions sufficient to form a bond between said GlcNAc and said EPO, wherein said GnT is selected from the group consisting of GnT I, GnT II, GnT III, GnT IV, GnT V and GnT VI.

54. The method of claim 53, wherein said mixture comprises one GnT
15 selected from the group consisting of GnT I, GnT II, GnT IV, GnT V and GnT VI.

55. The method of claim 54, wherein said GnT is GnT I.

56. The method of claim 54, wherein said GnT is GnT II.

57. The method of claim 39, wherein said glycoPEGylated EPO peptide
comprises at least one mono-antennary glycan.

58. The method of claim 39, wherein the sugar of said nucleotide sugar is
25 galactose and said glycosyltransferase is galactosyl transferase I (GalT I).

59. The method of claim 53, wherein before step (a) but after step (b):

(c) contacting said EPO peptide with a mixture comprising a nucleotide galactose
30 (Gal) and galactosyl transferase I (GalT I) under conditions sufficient to transfer galactose to said EPO peptide.

60. The method of claim 39, wherein in step (a), the sugar of said nucleotide sugar is sialic acid and said glycosyltransferase is a sialyltransferase.

5 61. The method of claim 60, wherein said sialic acid is N-acetylneuraminic acid (NAN).

62. The method of claim 60, wherein said sialyltransferase is selected from the group consisting of $\alpha(2,3)$ sialyltransferase, $\alpha(2,6)$ sialyltransferase and
10 (2,8)sialyltransferase.

63. A glycoPEGylated EPO peptide made by the method of claim 39.

64. A glycoPEGylated EPO peptide, said EPO peptide comprising the
15 sequence of SEQ ID NO:73.

65. A glycoPEGylated EPO peptide, said EPO peptide comprising the sequence of SEQ ID NO:73 and further comprising a mutation in said sequence.

20 66. A method of making a glycoPEGylated EPO peptide, said method comprising the steps of:

(a) contacting an EPO peptide with a mixture comprising a nucleotide sugar covalently linked to poly(ethylene glycol) and a glycosyltransferase under conditions sufficient to transfer said poly(ethylene glycol) to said EPO peptide, wherein said
25 glycosyltransferase is a fucosyltransferase.

67. The method of claim 66, wherein said fucosyltransferase is selected from the group consisting of fucosyltransferase I, fucosyltransferase III, fucosyltransferase IV, fucosyltransferase V, fucosyltransferase VI and fucosyltransferase VII.

30 68. A glycoPEGylated EPO peptide made by the method of claim 66.

69. The method of claim 66, wherein said EPO peptide is expressed in a CHO cell.

70. A method of treating a mammal having anemia, said method comprising administering to said mammal an EPO peptide having one or more glycans having a glycoconjugate molecule attached to said peptide, wherein said EPO peptide is administered in an amount effective to increase the hematocrit level in said mammal.

71. The method of claim 70, wherein said mammal is a human.

72. A method of providing erythropoietin therapy to a mammal, said method comprising administering an effective amount of a glycoPEGylated EPO peptide comprising an EPO peptide and at least one glycan and at least one poly(ethylene glycol) molecule covalently attached to said glycan, wherein said poly(ethylene glycol) molecule is added to said EPO peptide using a glycosyltransferase, wherein said EPO peptide is administered in an amount effective to increase the hematocrit level in said mammal.

73. The method of claim 72, wherein said mammal is a human.

74. A method of treating a mammal having anemia, said method comprising administering to said mammal a glycoPEGylated EPO peptide comprising an EPO peptide and at least one glycan and at least one poly(ethylene glycol) molecule covalently attached to said glycan, wherein said poly(ethylene glycol) molecule is added to said EPO peptide using a glycosyltransferase, wherein said EPO peptide is administered in an amount effective to increase the hematocrit level in said mammal..

75. The method of claim 74, wherein said mammal is a human.

76. The method of claim 75, wherein said anemia is associated with chemotherapy.

77. A method of treating a kidney dialysis patient, said method comprising
5 administering to said patient a glycoPEGylated EPO peptide comprising an EPO peptide and
at least one glycan and at least one poly(ethylene glycol) molecule covalently attached to said
glycan, wherein said poly(ethylene glycol) molecule is added to said EPO peptide using a
glycosyltransferase, wherein said EPO peptide is administered in an amount effective to
increase the hematocrit level in said patient.

10

1/498

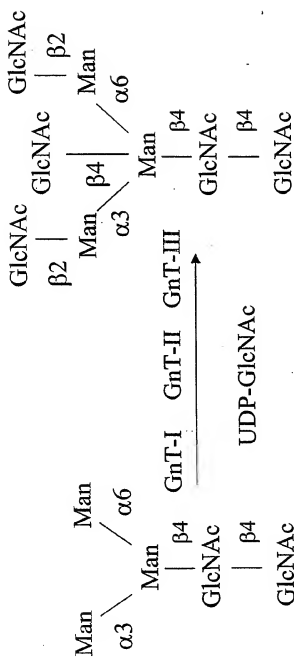
Trimannosyl core with
Bisecting GlcNAc

FIG. 1

Trimannosyl core

2/498


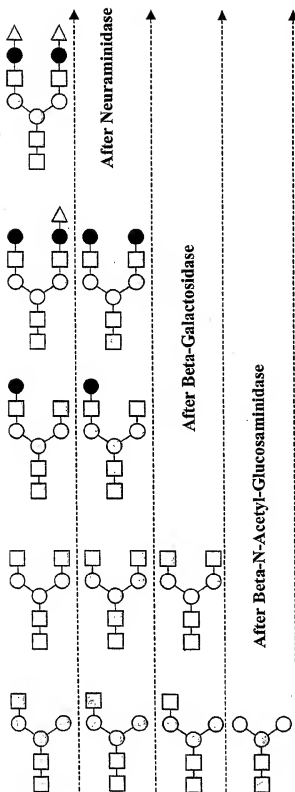


FIG. 2

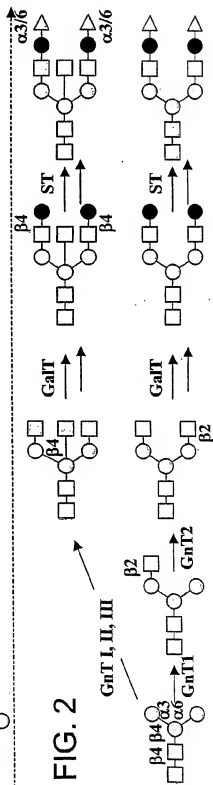

3/498

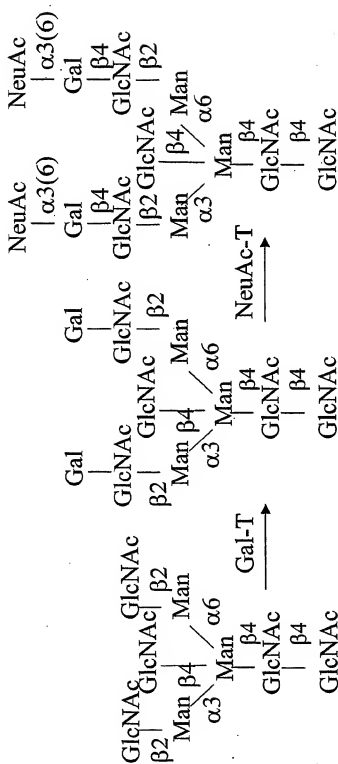


FIG. 3

4/498

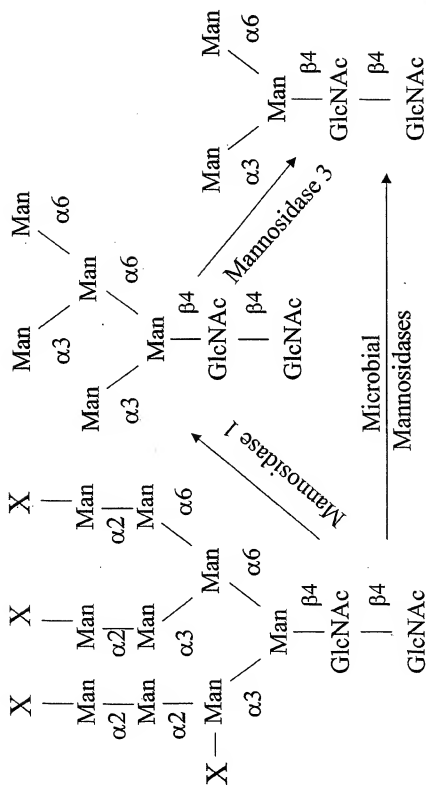


FIG. 4

5/498

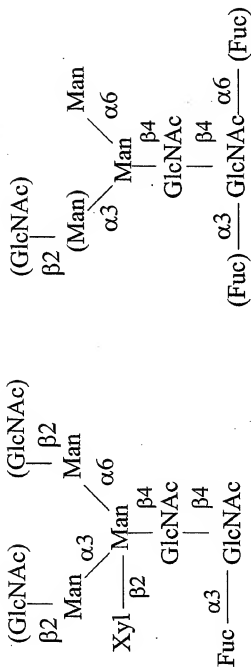


FIG. 5

FIG. 6

6/498

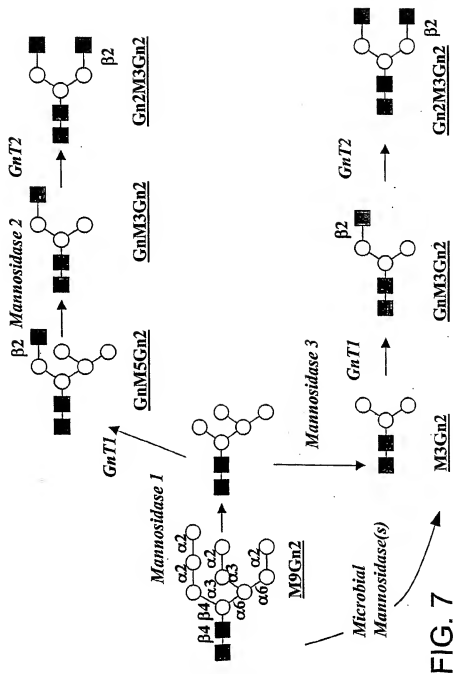


FIG. 7

7/498

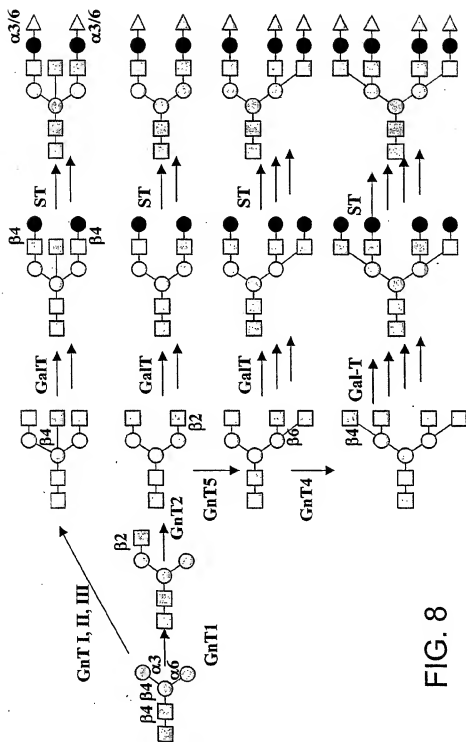


FIG. 8

8/498

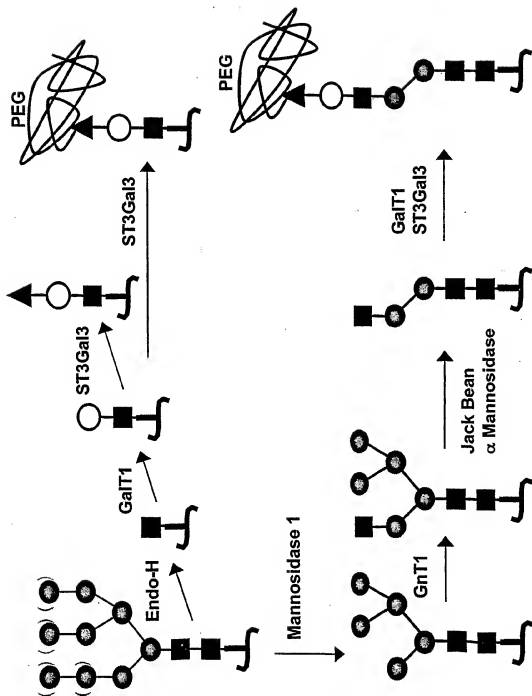


FIG. 9

9/498

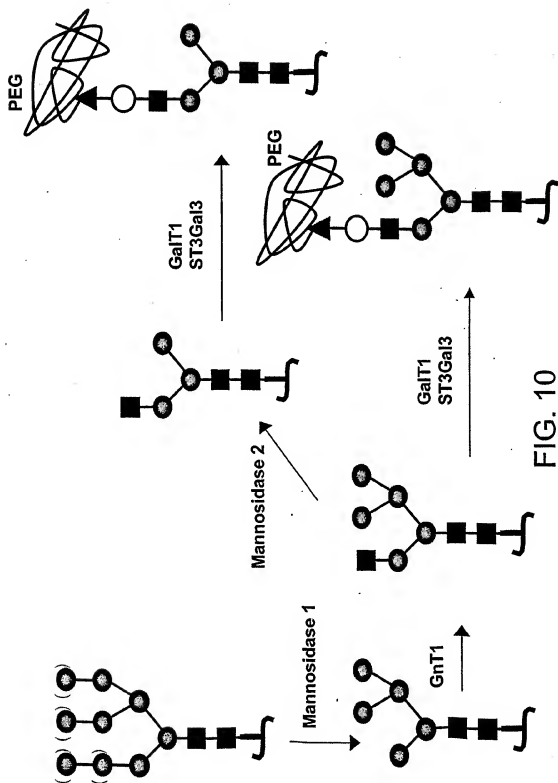


FIG. 10

10/498

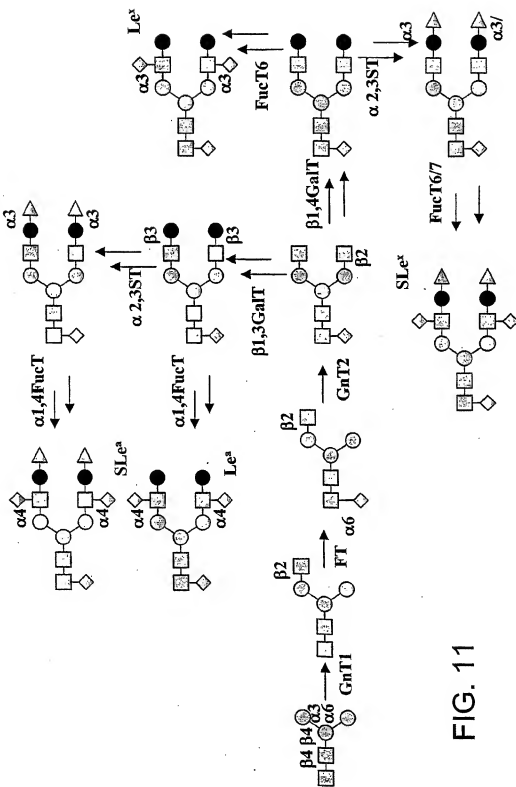


FIG. 11

11/498

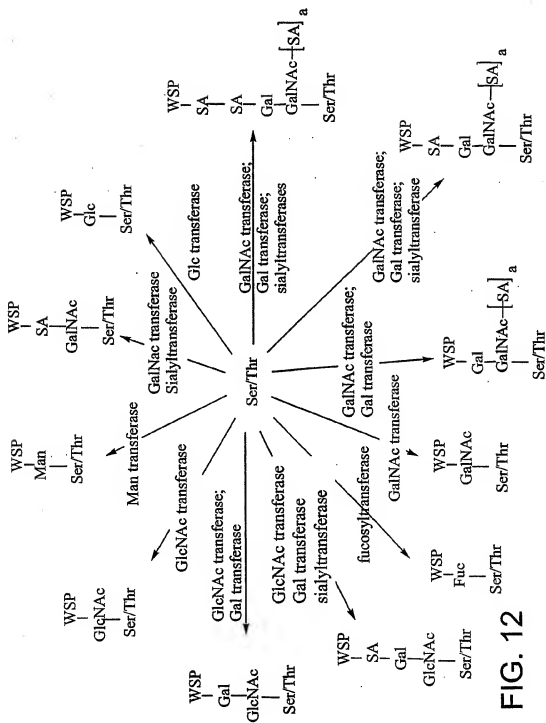


FIG. 12

12/498

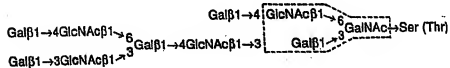
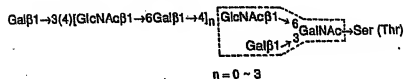
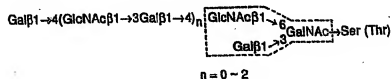
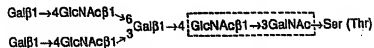
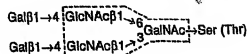
Core 1**Core 2****Core 3****Core 4**

FIG. 13

14/498

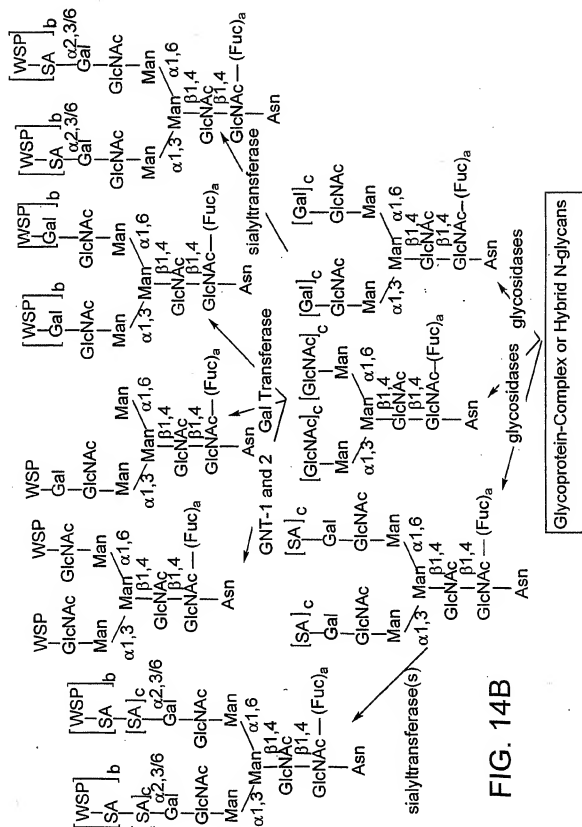


FIG. 14B

17/498

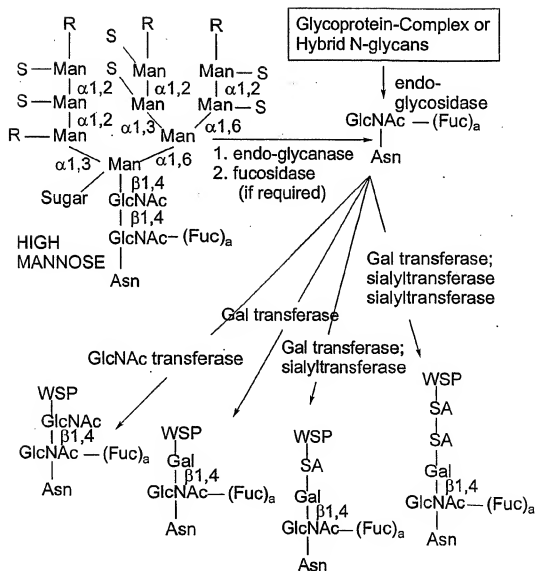


FIG. 17

22/498

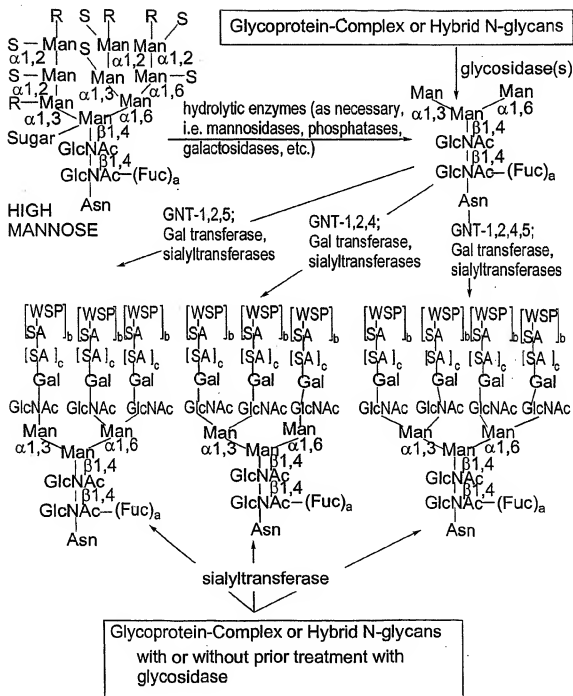


FIG. 20

25/498

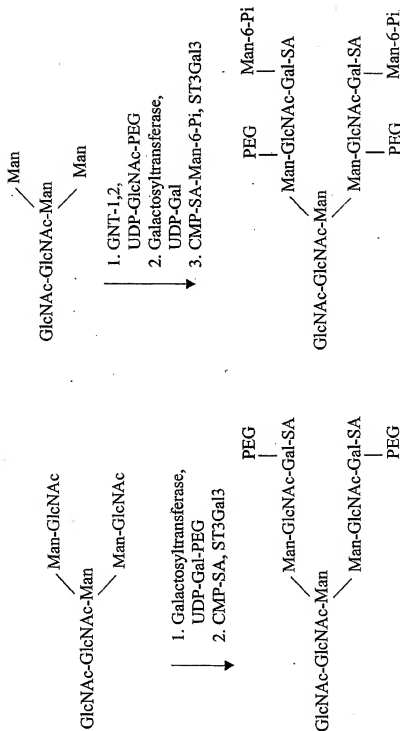


FIG. 23A

FIG. 23B

26/498

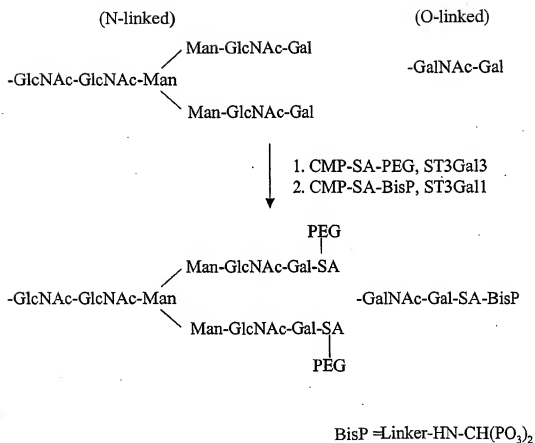


FIG. 23C

27/498

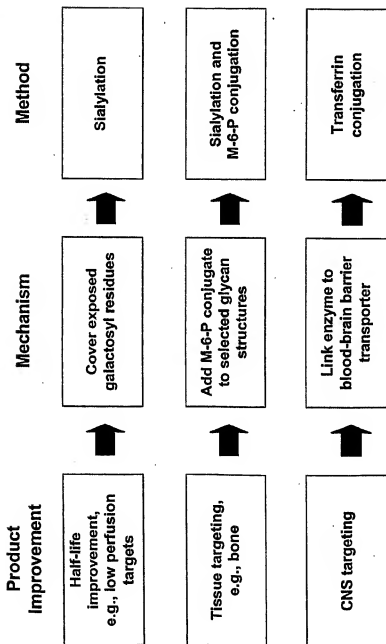


FIG. 24

28/498

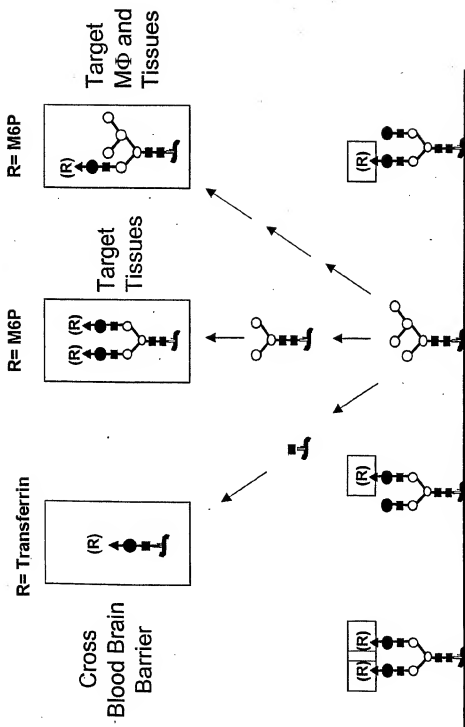


FIG. 25

29/498

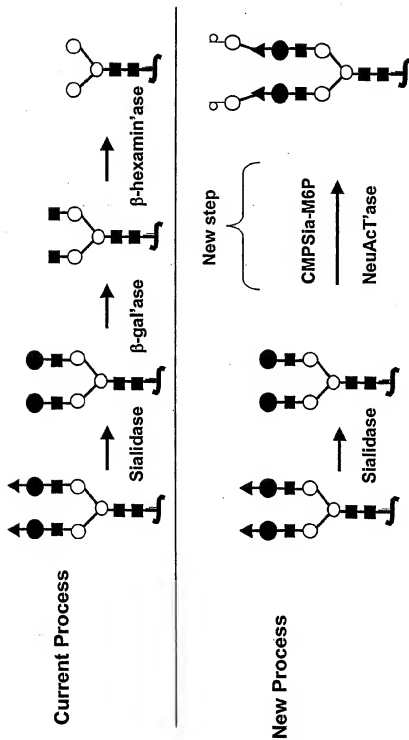


FIG. 26

30/498

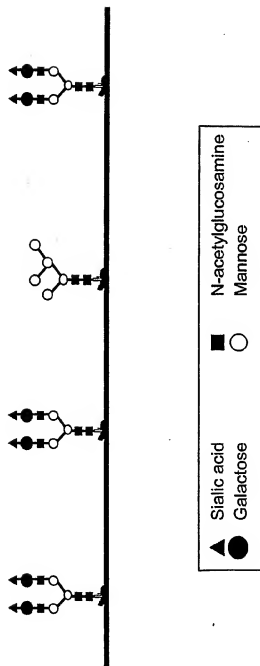


FIG. 27

31/498

12AP1/E5 -- Viventia Biotech
 1964 -- Aventis
 20K growth hormone -- AMUR
 28P6/E6 -- Viventia Biotech
 3-Hydroxyphthaloyl-beta-lactoglobulin --
 4-IBB ligand gene therapy --
 64-Cu MAb conjugate TETA-1A3 --
 Mallinckrodt Institute of Radiology
 64-Cu MAb conjugate TETA-cT84.66
 64-Cu Trastuzumab TETA conjugate --
 Genentech
 A 200 -- Amgen
 A10255 -- Eli Lilly
 A1PDX -- Hedral Therapeutics
 A6 -- Angstrom
 aaAT-III -- Genzyme
 Abciximab -- Centocor
 ABI.001 -- Atlantic BioPharmaceuticals
 ABT-828 -- Abbott
 Accutin
 Actinohivin
 activin -- Biotech Australia, Human
 Therapeutics, Curis
 AD 439 -- Tanox
 AD 519 -- Tanox
 Adalimumab -- Cambridge Antibody Tech.
 Adenocarcinoma vaccine -- Biomira -- NIS
 Adenosine deaminase -- Enzond
 Adenosine A2B receptor antagonists --
 Adenosine Therapeutics
 ADP-001 -- Axis Genetics
 AF 13948 -- Affymax
 Afelimomab -- Knoll
 AFP-SCAN -- Immunomedics
 AG 2195 -- Corixa
 agalsidase alfa -- Transkaryotic Therapies
 agalsidase beta -- Genzyme
 AGENT-- Antisoma
 AI 300 -- Autolmmune
 AI-101 -- Teva
 AI-102 -- Teva
 AI-201 -- Autolmmune
 AI-301 -- Autolmmune
 AIDS vaccine -- ANRS, CIBG, Hese
 Biomed, Hollis-Eden, Rome, United
 Biomedical, American Home Products,
 Maxygen
 airway receptor ligand -- IC Innovations
 AJW 2 -- Ajinomoto
 AK 30 NGF -- Alkermes
 Albuferon -- Human Genome Sciences
 albumin -- Biogen, DSM Anti-Infectives,
 Genzyme Transgenics, PPL Therapeutics,
 TranXenoGen, Welfide Corp.
 aldesleukin -- Chiron
 alefacept -- Biogen
 Alemtuzumab
 Allergy therapy -- ALK-Abello/Maxygen,
 ALK-Abello/RP Scherer
 allergy vaccines -- Allergy Therapeutics
 Alnidofibatide -- Aventis Pasteur
 Alnorine -- SRC VB VECTOR
 ALP 242 -- Gruenenthal
 Alpha antitrypsin -- Arriva/Hyland
 Immuno/ProMetic/Protease Sciences
 Alpha-1 antitrypsin -- Cutter, Bayer, PPL
 Therapeutics, Profile, ZymoGenetics,
 Arriva
 Alpha-1 protease inhibitor -- Genzyme
 Transgenics, Welfide Corp.
 Alpha-galactose fusion protein --
 Immunomedics
 Alpha-galactosidase A -- Research
 Corporation Technologies, Genzyme
 Alpha-glucosidase -- Genzyme, Novazyme
 Alpha-lactalbumin
 Alpha-L-iduronidase -- Transkaryotic
 Therapies, BioMarin
 alteplase -- Genentech
 alvircept sudotox -- NIH
 ALX-0600, a GLP-2 agonist -- NPS Allelix
 Corp.

FIG. 28A

32/498

ALX1-11 --sNPS Pharmaceuticals
 Alzheimer's disease gene therapy
 AM-133 -- AMRAD
 Amb a 1 immunostim conj. -- Dynavax
 AMD 3100 -- AnorMED -- NIS
 AMD 3465 -- AnorMED -- NIS
 AMD 3465 -- AnorMED -- NIS
 AMD Fab -- Genentech
 Amediplase -- Menarini, Novartis
 AM-F9
 Amoebiasis vaccine
 Amphiregulin -- Octagene
 anakinra -- Amgen
 analgesic -- Nobex
 aneastim -- Amgen
 Anergix.RA -- Corixa, Organon
 Angiocidin -- InKine
 angiogenesis inhibitors -- ILEX
 AngioMab -- Antisoma
 Angiopoietins -- Regeneron/Procter &
 Gamble
 angiostatin -- EntreMed
 Angiostatin/endostatin gene therapy --
 Genetix Pharmaceuticals
 angiotensin-II, topical -- Maret
 Anthrax -- EluSys Therapeutics/US Army
 Medical Research Institute
 Anthrax vaccine
 Anti platelet-derived growth factor D human
 monoclonal antibodies -- CuraGen
 Anti-17-1A Mab 3622W94 --
 GlaxoSmithKline
 Anti-2C4 Mab -- Genentech
 anti-4-1BB monoclonal antibodies -- Bristol-
 Myers Squibb
 Anti-Adhesion Platform Tech. -- CytoVax
 Anti-adipocyte Mab -- Cambridge Antibody
 Tech./ObeSys
 antiallergics -- Maxygen
 antiallergy vaccine -- Acambis
 Anti-alpha-4-integrin Mab
 Anti-alphavβ3 integrin Mab -- Applied
 Molecular Evolution
 Anti-angiogenesis monoclonal antibodies --
 KS Biomedix/Schering AG
 Anti-B4 Mab-DC1 conjugate -- ImmunoGen
 Anti-B7 antibody PRIMATIZED -- IDEC
 Anti-B7-1 Mab 16-10A1
 Anti-B7-1 Mab 1G10
 Anti-B7-2 Mab GL-1
 Anti-B7-2-gelonin immunotoxin --
 Antibacterials/antifungals --
 Diversa/IntraBiotics
 Anti-beta-amyloid monoclonal antibodies --
 Cambridge Antibody Tech., Wyeth-Ayerst
 Anti-BLyS antibodies -- Cambridge
 Antibody Tech./Human Genome Sciences
 Antibody-drug conjugates -- Seattle
 Genetics/Eos
 Anti-C5 Mab BB5-1 -- Alexion
 Anti-C5 Mab N19-8 -- Alexion
 Anti-C8 Mab
 anticancer cytokines -- BioPulse
 anticancer matrix -- Telios Integra
 Anticancer monoclonal antibodies -- ARIUS,
 Immunex
 anticancer peptides -- Maxygen, Micrologix
 Anticancer prodrug Tech. -- Alexion
 Antibody Technologies
 anticancer Troy-Bodies -- Affite -- Affitech
 anticancer vaccine -- NIH
 anticancers -- Epimmune
 Anti-CCR5/CXCR4 sheep Mab -- KS
 Biomedix Holdings
 Anti-CD11a Mab KBA --
 Anti-CD11a Mab M17
 Anti-CD11a Mab TA-3 --
 Anti-CD11a Mab WT.1 --
 Anti-CD11b Mab -- Pharmacia
 Anti-CD11b Mab LM2
 Anti-CD154 Mab -- Biogen
 Anti-CD16-anti-CD30 Mab -- Biotest

FIG. 28B

33/498

Anti-CD18 MAb -- Pharmacia	Anti-CD4 MAb -- Centocor, IDEC
Anti-CD19 MAb B43 --	Pharmaceuticals, Xenova Group
Anti-CD19 MAb -liposomal sodium butyrate conjugate --	Anti-CD4 MAb 16H5
Anti-CD147	Anti-CD4 MAb 4162W94 -- GlaxoSmithKline
Anti-CD19 MAb-saporin conjugate --	Anti-CD4 MAb B-F5 -- Diaclone
Anti-CD19-dsFv-PE38-immunotoxin --	Anti-CD4 MAb GK1-5
Anti-CD2 MAb 12-15 --	Anti-CD4 MAb KT6
Anti-CD2 MAb B-E2 -- Diaclone	Anti-CD4 MAb OX38
Anti-CD2 MAb OX34 --	Anti-CD4 MAb PAP conjugate -- Bristol-Myers Squibb
Anti-CD2 MAb OX54 --	Anti-CD4 MAb RIB 5-2
Anti-CD2 MAb OX55 --	Anti-CD4 MAb W3/25
Anti-CD2 MAb RM2-1	Anti-CD4 MAb YTA 3.1.2
Anti-CD2 MAb RM2-2	Anti-CD4 MAb YTS 177-9
Anti-CD2 MAb RM2-4	Anti-CD40 ligand MAb 5c8 -- Biogen
Anti-CD20 MAb BCA B20	Anti-CD40 MAb
Anti-CD20-anti-Fc alpha RI bispecific MAb -- Medarex, Tenovus	Anti-CD40 MAb 5D12 -- Tanox
Anti-CD22 MAb-saporin-6 complex --	Anti-CD44 MAb A3D8
Anti-CD3 Immunotoxin --	Anti-CD44 MAb GKWA3
Anti-CD3 MAb 145-2C11 -- Pharming	Anti-CD44 MAb IM7
Anti-CD3 MAb CD4IgG conjugate -- Genentech	Anti-CD44 MAb KM81
Anti-CD3 MAb humanised -- Protein Design, RW Johnson	Anti-CD44 variant monoclonal antibodies -- Corixa/Hebrew University
Anti-CD3 MAb WT32	Anti-CD45 MAb BC8-I-131
Anti-CD3 MAb-ricin-chain-A conjugate --	Anti-CD45RB MAb
Anti-CD3 MAb-xanthine-oxidase conjugate --	Anti-CD48 MAb HuLy-m3
Anti-CD30 MAb BerH2 -- Medac	Anti-CD48 MAb WM-63
Anti-CD30 MAb-saporin conjugate	Anti-CD5 MAb -- Becton Dickinson
Anti-CD30-scFv-ETA'-immunotoxin	Anti-CD5 MAb OX19
Anti-CD38 MAb AT13/5	Anti-CD6 MAb
Anti-CD38 MAb-saporin conjugate	Anti-CD7 MAb-PAP conjugate
Anti-CD3-anti-CD19 bispecific MAb	Anti-CD7 MAb-ricin-chain-A conjugate
Anti-CD3-anti-EGFR MAb	Anti-CD8 MAb -- Amerimmune, Cytodyn, Becton Dickinson
Anti-CD3-anti-interleukin-2-receptor MAb	Anti-CD8 MAb 2-43
Anti-CD3-anti-MOV18 MAb -- Centocor	Anti-CD8 MAb OX8
Anti-CD3-anti-SCLC bispecific MAb	Anti-CD80 MAb P16C10 -- IDEC
Anti-CD4 idiotype vaccine	Anti-CD80 MAb P7C10 -- ID Vaccine
	Anti-CD8-idarubicin conjugate
	Anti-CEA MAb CE-25
	Anti-CEA MAb MN 14 -- Immunomedics

FIG. 28C

34/498

- Anti-CEA MAb MN14-PE40 conjugate -- Immunomedics
- Anti-CEA MAb T84.66-interleukin-2 conjugate
- Anti-CEA sheep MAb -- KS Biomedix Holdings
- Anti-cell surface monoclonal antibodies -- Cambridge Antibody Tech. /Pharmacia
- Anti-c-erbB2-anti-CD3 bifunctional MAb -- Otsuka
- Anti-CMV MAb -- Scotgen
- Anti-complement
- Anti-CTLA-4 MAb
- Anti-EGFR catalytic antibody -- Hersed Biomed
- anti-EGFR immunotoxin -- IVAX
- Anti-EGFR MAb -- Abgenix
- Anti-EGFR MAb 528
- Anti-EGFR MAb KSB 107 -- KS Biomedix
- Anti-EGFR MAb-DM1 conjugate -- ImmunoGen
- Anti-EGFR MAb-LA1 --
- Anti-EGFR sheep MAb -- KS Biomedix
- Anti-FAP MAb F19-I-131
- Anti-Fas IgM MAb CH11
- Anti-Fas MAb Jo2
- Anti-Fas MAb RK-8
- Anti-Fit-1 monoclonal antibodies -- ImClone
- Anti-fungal peptides -- State University of New York
- antifungal tripeptides -- BTG
- Anti-ganglioside GD2 antibody-interleukin-2 fusion protein -- Lexigen
- Anti-GM2 MAb -- Kyowa
- Anti-GM-CSF receptor monoclonal antibodies -- AMRAD
- Anti-gp130 MAb -- Tosoh
- Anti-HCA monoclonal antibodies -- AltaRex/Epigen
- Anti-hCG antibodies -- Abgenix/AVI BioPharma
- Anti-heparanase human monoclonal antibodies -- Oxford Glycosciences/Medarex
- Anti-hepatitis C virus human monoclonal antibodies -- XTL Biopharmaceuticals
- Anti-HER-2 antibody gene therapy
- Anti-herpes antibody -- Epicyte
- Anti-HIV antibody -- Epicyte
- anti-HIV catalytic antibody -- Hersed Biomed
- anti-HIV fusion protein -- Idun
- anti-HIV proteins -- Cangene
- Anti-HM1-24 MAb -- Chugai
- Anti-hrR3 MAb
- Anti-Human-Carcinoma-Antigen MAb -- Epicyte
- Anti-ICAM-1 MAb -- Boehringer Ingelheim
- Anti-ICAM-1 MAb 1A-29 -- Pharmacia
- Anti-ICAM-1 MAb HA58
- Anti-ICAM-1 MAb YN1/1.7.4
- Anti-ICAM-3 MAb ICM3 -- ICOS
- Anti-idiotypic breast cancer vaccine 11D10
- Anti-idiotypic breast cancer vaccine ACA14C5 --
- Anti-idiotypic cancer vaccine -- ImClone Systems/Merck KGaA ImClone, Viventia Biotech
- Anti-idiotypic cancer vaccine 1A7 -- Titan
- Anti-idiotypic cancer vaccine 3H1 -- Titan
- Anti-idiotypic cancer vaccine TriAb -- Titan
- Anti-idiotypic Chlamydia trachomatis vaccine
- Anti-idiotypic colorectal cancer vaccine -- Novartis
- Anti-idiotypic colorectal cancer vaccine -- Onyvox
- Anti-idiotypic melanoma vaccine -- IDEC Pharmaceuticals
- Anti-idiotypic ovarian cancer vaccine ACA 125
- Anti-idiotypic ovarian cancer vaccine AR54 - AltaRex

FIG. 28D

35/498

Anti-idiotypic ovarian cancer vaccine CA-125 – AltaRex, Biomira	Anti-L-selectin monoclonal antibodies -- Protein Design Labs, Abgenix, Stanford University
Anti-IgE catalytic antibody -- Hesed Biomed	Anti-MBL monoclonal antibodies -- Alexion/Brigham and Women's Hospital
Anti-IgE MAb E26 -- Genentech	Anti-MHC monoclonal antibodies
Anti-IGF-1 MAb	Anti-MIF antibody humanised -- IDEC, Cytokine PharmaSciences
anti-inflammatory -- GeneMax	Anti-MRSA/VRSA sheep MAb -- KS Biomedix Holdings
anti-inflammatory peptide -- BTG	Anti-mu MAb -- Novartis
anti-integrin peptides -- Burnha	Anti-MUC-1 MAb
Anti-interferon-alpha-receptor MAb 64G12 -- Pharma Pacific Management	Anti-MUC 18
Anti-interferon-gamma MAb -- Protein Design Labs	Anti-Nogo-A MAb IN1
Anti-interferon-gamma polyclonal antibody - Advanced Biotherapy	Anti-nuclear autoantibodies -- Procyon
Anti-interleukin-10 MAb --	Anti-ovarian cancer monoclonal antibodies - Dompe
Anti-interleukin-12 MAb --	Anti-p185 monoclonal antibodies
Anti-interleukin-1-beta polyclonal antibody -- R&D Systems	Anti-p43 MAb
Anti-interleukin-2 receptor MAb 2A3	Antiparasitic vaccines
Anti-interleukin-2 receptor MAb 33B3-1 -- Immunotech	Anti-PDGF/bFGF sheep MAb -- KS Biomedix
Anti-interleukin-2 receptor MAb ART-18	Anti-properdin monoclonal antibodies -- Abgenix/Gliatech
Anti-interleukin-2 receptor MAb LO-Tact-1	Anti-PSMA (prostate specific membrane antigen)
Anti-interleukin-2 receptor MAb Mikbeta1	Anti-PSMA MAb J591 -- BZL Biologics
Anti-interleukin-2 receptor MAb NDS61	Anti-Rev MAb gene therapy --
Anti-interleukin-4 MAb 11B11	Anti-RSV antibodies -- Epicyte, Intracell
Anti-interleukin-5 MAb -- Wallace Laboratories	Anti-RSV monoclonal antibodies -- Medarex/MedImmune, Applied Molecular Evolution/MedImmune
Anti-interleukin-6 MAb -- Centocor, Diacorne, Pharmadigm	Anti-RSV MAb, inhalation -- Alkermes/MedImmune
Anti-interleukin-8 MAb -- Abgenix	Anti-RT gene therapy
Anti-interleukin-8 MAb -- Xenotech	Antisense K-ras RNA gene therapy
Anti-JL1 MAb	Anti-SF-25 MAb
Anti-Klebsiella sheep MAb -- KS Biomedix Holdings	Anti-sperm antibody -- Epicyte
Anti-Laminin receptor MAb-liposomal doxorubicin conjugate	Anti-Tac(Fv)-PE38 conjugate
Anti-LCG MAb -- Cytoclonal	Anti-TAPA/CD81 MAb AMP1
Anti-lipopolysaccharide MAb -- VitaResc	Anti-tat gene therapy

36/498

Anti-TCR-alphabeta MAb H57-597
 Anti-TCR-alphabeta MAb R73
 Anti-tenascin MAb BC-4-I-131
 Anti-TGF-beta human monoclonal antibodies -- Cambridge Antibody Tech., Genzyme
 Anti-TGF-beta MAb 2G7 -- Genentech
 Antithrombin III -- Genzyme Transgenics, Aventis, Bayer, Behringwerke, CSL, Myriad
 Anti-Thy1 MAb
 Anti-Thy1.1 MAb
 Anti-tissue factor/factor VIIA sheep MAb -- KS Biomedix
 Anti-TNF monoclonal antibodies -- Centocor, Chiron, Peptech, Pharacia, Sero
 Anti-TNF sheep MAb -- KS Biomedix Holdings
 Anti-TNFalpha MAb -- Genzyme
 Anti-TNFalpha MAb B-C7 -- Diaclone
 Anti-tooth decay MAb -- Planet BioTech.
 Anti-TRAIL receptor-1 MAb -- Takeda
 Antitumour RNases -- NIH
 Anti-VCAM MAb 2A2 -- Alexion
 Anti-VCAM MAb 3F4 -- Alexion
 Anti-VCAM-1 MAb
 Anti-VEC MAb -- ImClone
 Anti-VEGF MAb -- Genentech
 Anti-VEGF MAb 2C3
 Anti-VEGF sheep MAb -- KS Biomedix Holdings
 Anti-VLA-4 MAb HP1/2 -- Biogen
 Anti-VLA-4 MAb PS/2
 Anti-VLA-4 MAb R1-2
 Anti-VLA-4 MAb TA-2
 Anti-VAP-1 human MAb
 Anti-VRE sheep MAb -- KS Biomedix Holdings
 ANUP -- TranXenoGen
 ANUP-1 -- Pharis
 AOP-RANTES -- Senetek
 Apan-CH -- Praecis Pharmaceuticals
 APC-8024 -- Demegen
 ApoA-1 -- Milano, Pharmacia
 Apogen -- Alexion
 apolipoprotein A1 -- Avanir
 Apolipoprotein E -- Bio-Tech. General
 Applaggin -- Biogen
 aprotinin -- ProdiGene
 APT-070C -- AdProTech
 AR 177 -- Aronex Pharmaceuticals
 AR 209 -- Aronex Pharmaceuticals, Antigenics
 AR545C
 ARGENT gene delivery systems -- ARIAD
 Arresten
 ART-123 -- Asahi Kasei
 arylsulfatase B -- BioMarin
 Arylsulfatase B, Recombinant human -- BioMarin
 AS 1051 -- Ajinomoto
 ASI-BCL -- Intracell
 Asparaginase - Merck
 ATL-101 -- Alizyme
 Atrial natriuretic peptide -- Pharis
 Aurintricarboxylic acid-high molecular weight
 Autoimmune disorders -- GPC
 Biotech/MorphoSys
 Autoimmune disorders and transplant rejection -- Bristol-Myers Squibb/Genzyme
 Tra
 Autoimmune disorders/cancer -- Abgenix/Chiron, CuraGen
 Autotaxin
 Avicidin -- NeoRx
 axogenesis factor-1 -- Boston Life Sciences
 Axokine -- Regeneron
 B cell lymphoma vaccine -- Biomira
 B7-1 gene therapy --
 BABS proteins -- Chiron

FIG. 28F

37/498

BAM-002 -- Novelos Therapeutics
 Basiliximab (anti CD25 MAb) -- Novartis
 Bay-16-9996 -- Bayer
 Bay-39-9437 -- Bayer
 Bay-50-4798 -- Bayer
 BB-10153 -- British Biotech
 BBT-001 -- Bolder BioTech.
 BBT-002 -- Bolder BioTech.
 BBT-003 -- Bolder BioTech.
 BBT-004 -- Bolder BioTech.
 BBT-005 -- Bolder BioTech.
 BBT-006 -- Bolder BioTech.
 BBT-007 -- Bolder BioTech.
 BCH-2763 -- Shire
 BCSF -- Millenium Biologix
 BDNF -- Regeneron -- Amgen
 Becapiermin -- Johnson & Johnson, Chiron
 Bectumomab -- Immunomedics
 Beriplast -- Aventis
 Beta-adrenergic receptor gene therapy --
 University of Arkansas
 bFGF -- Scios
 BI 51013 -- Behringwerke AG
 BIBH 1 -- Boehringer Ingelheim
 BIM-23190 -- Beaufour-Ipsen
 birch pollen immunotherapy -- Pharmacia
 bispecific fusion proteins -- NIH
 Bispecific MAb 2B1 -- Chiron
 Bitistatin
 BIWA 4 -- Boehringer Ingelheim
 blood substitute -- Northfield, Baxter Intl.
 BLP-25 -- Biomira
 BLS-0597 -- Boston Life Sciences
 BLyS -- Human Genome Sciences
 BLyS radiolabelled -- Human Genome
 Sciences
 BM 06021 -- Boehringer Mannheim
 BM-202 -- BioMarin
 BM-301 -- BioMarin
 BM-301 -- BioMarin
 BM-302 -- BioMarin
 BMP 2 -- Genetics Institute/Medtronic-
 Sofamor Danek, Genetics Institute/
 Collagenesis, Genetics
 Institute/Yamanouch
 BMP 2 gene therapy
 BMP 52 -- Aventis Pasteur, Biopharm
 BMP-2 -- Genetics Institute
 BMS 182248 -- Bristol-Myers Squibb
 BMS 202448 -- Bristol-Myers Squibb
 bone growth factors -- IsoTis
 BPC-15 -- Pfizer
 brain natriuretic peptide --
 Breast cancer -- Oxford
 GlycoSciences/Medarex
 Breast cancer vaccine -- Therion Biologics,
 Oregon
 BSSL -- PPL Therapeutics
 BST-2001 -- BioStratum
 BST-3002 -- BioStratum
 BTI 322 --
 butyrylcholinesterase -- Shire
 C 6822 -- COR Therapeutics
 C1 esterase inhibitor -- Pharming
 C3d adjuvant -- AdProTech
 CAB-2.1 -- Millennium
 calcitonin -- Inhale Therapeutics Systems,
 Aventis, Genetronics, TranXenoGen,
 Unigene, Rhone Poulenc Rohrer
 calcitonin -- oral -- Nobex, Emisphere,
 Pharmaceutical Discovery
 Calcitonin gene-related peptide -- Asahi
 Kasei -- Unigene
 calcitonin, human -- Suntory
 calcitonin, nasal -- Novartis, Unigene
 calcitonin, Panoderm -- Elan
 calcitonin, Peptitrol -- Shire
 calcitonin, salmon -- Therapicon
 calin -- Biopharm
 Calphobindin I
 calphobindin I -- Kowa
 calreticulin -- NYU

FIG. 28G

38/498

Campath-1G
 Campath-1M
 cancer therapy -- Cangene
 cancer vaccine -- Aixlie, Aventis Pasteur,
 Center of Molecular Immunology, YM
 BioSciences, Cytos, Genzyme,
 Transgenics, GlobelImmune, Igeneon,
 ImClone, Virogenetics, InterCell, Iomai,
 Jenner Biotherapies, Memorial Sloan-
 Kettering Cancer Center, Sydney Kimmel
 Cancer Center, Novavax, Protein
 Sciences, Argonex, SIGA
 Cancer vaccine ALVAC-CEA B7.1 --
 Aventis Pasteur/Therion Biologics
 Cancer vaccine CEA-TRICOM -- Aventis
 Pasteur/Therion Biologics
 Cancer vaccine gene therapy -- Cantab
 Pharmaceuticals
 Cancer vaccine HER-2/neu -- Corixa
 Cancer vaccine THERATOPE -- Biomira
 cancer vaccine, PolyMASC -- Valentis
 Candida vaccine -- Corixa, Inhibitex
 Canstatin -- ILEX
 CAP-18 -- Panorama
 Cardiovascular gene therapy -- Collateral
 Therapeutics
 carperitide -- Suntory
 Casocidin-1 -- Pharis
 CAT 152 -- Cambridge Antibody Tech.
 CAT 192 -- Cambridge Antibody Tech.
 CAT 213 -- Cambridge Antibody Tech.
 Catalase-- Enzon
 Cat-PAD -- Circassia
 CB 0006 -- Celltech
 CCK(27-32)-- Akzo Nobel
 CCR2-64I -- NIH
 CD, Procept -- Paligent
 CD154 gene therapy
 CD39 -- Immunex
 CD39-L2 -- Hyseq
 CD39-L4 -- Hyseq
 CD4 fusion toxin -- Senetek
 CD4 IgG -- Genentech
 CD4 receptor antagonists --
 Pharmacocepeia/Progenics
 CD4 soluble -- Progenics
 CD4, soluble -- Genzyme Transgenics
 CD40 ligand -- Immunex
 CD4-ricin chain A -- Genentech
 CD59 gene therapy -- Alexion
 CD8 TIL cell therapy -- Aventis Pasteur
 CD8, soluble -- Avidex
 CD95 ligand -- Roche
 CDP 571 -- Celltech
 CDP 850 -- Celltech
 CDP-860 (PEG-PDGF MAb) -- Celltech
 CDP 870 -- Celltech
 CDS-1 -- Ernest Orlando
 Cedelizumab -- Ortho-McNeil
 Cetermin -- Insmed
 CETP vaccine -- Avant
 Cetorelix
 Cetuximab
 CGH 400 -- Novartis
 CGP 42934 -- Novartis
 CGP 51901 -- Tanox
 CGRP -- Unigene
 CGS 27913 -- Novartis
 CGS 32359 -- Novartis
 Chagas disease vaccine -- Corixa
 chemokines -- Immune Response
 CHH 380 -- Novartis
 chitinase -- Genzyme, ICOS
 Chlamydia pneumoniae vaccine -- Antex
 Biologics
 Chlamydia trachomatis vaccine -- Antex
 Biologics
 Chlamydia vaccine -- GlaxoSmithKline
 Cholera vaccine CVD 103-HgR -- Swiss
 Serum and Vaccine Institute Berne
 Cholera vaccine CVD 112 -- Swiss Serum
 and Vaccine Institute Berne

FIG. 28H

39/498

Cholera vaccine inactivated oral -- SBL	CRL 1605 -- CytRx
Vaccin	CS-560 -- Sankyo
Chrysalin -- Chrysalis BioTech.	CSF -- ZymoGenetics
CI-782 -- Hitachi Kase	CSF-G -- Hangzhou, Dong-A, Hanmi
Ciliary neurotrophic factor -- Fidia, Roche	CSF-GM -- Cangene, Hunan, LG Chem
CIM project -- Active Biotech	CSF-M -- Zarix
CL 329753 -- Wyeth-Ayerst	CT 1579 -- Merck Frosst
CL22, Cobra -- ML Laboratories	CT 1786 -- Merck Frosst
Clenoliximab -- IDEC	CT-112 ^A -- BTG
Clostridium difficile antibodies -- Epicycle	CTB-134L -- Xenova
clotting factors -- Octagene	CTC-111 -- Kaketsuken
CMB 401 -- Celltech	CTGF -- FibroGen
CNTF -- Sigma-Tau	CTLA4-Ig -- Bristol-Myers Squibb
Cocaine abuse vaccine -- Cantab,	CTLA4-Ig gene therapy --
ImmuLogic, Scripps	CTP-37 -- AVI BioPharma
coccidiomycosis vaccine -- Arizo	C-type natriuretic peptide -- Suntory
collagen -- Type I -- Pharming	CVS 995 -- Corvas Intl.
Collagen formation inhibitors -- FibroGen	CX 397 -- Nikko Kyodo
Collagen/hydroxyapatite/bone growth factor	CY 1747 -- Epimmune
-- Aventis Pasteur, Biopharm, Orquest	CY 1748 -- Epimmune
collagenase -- BioSpecifics	Cyanovirin-N
Colorectal cancer vaccine -- Wistar Institute	Cystic fibrosis therapy -- CBR/IVAX
Component B, Recombinant -- Sero	CYT 351
Connective tissue growth factor inhibitors --	cytokine Traps -- Regeneron
FibroGen/Taisho	cytokines -- Enzon, Cytodonal
Contortrostatin	Cytomegalovirus glycoprotein vaccine --
contraceptive vaccine -- Zonagen	Chiron, Aquila Biopharmaceuticals,
Contraceptive vaccine hCG	Aventis Pasteur, Virogenetics
Contraceptive vaccine male reversible --	Cytomegalovirus vaccine live -- Aventis
IMMUCON	Pasteur
Contraceptive vaccine zona pellucida --	Cytosine deaminase gene therapy --
Zonagen	GlaxoSmithKline
Copper-64 labelled Mab TETA-1A3 -- NCI	DA-3003 -- Dong-A
Coralyne	DAB389interleukin-6 -- Senetek
Corsevin M	DAB389interleukin-7
C-peptide analogues -- Schwarz	DAC:GLP-2 -- ConjuChem, Inc.
CPI-1500 -- Consensus	Daclizumab (anti-IL2R Mab) -- Protein
CRF -- Neurobiological Tech.	Design Labs
cRGDfV pentapeptide --	DAMP ^A -- Incyte Genomics
CRL 1095 -- CytRx	Daniplestim -- Pharmacia
CRL 1336 -- CytRx	darbepoetin alfa -- Amgen

40/498

DBI-3019 -- Diabetogen
 DCC -- Genzyme
 DDF -- Hyseq
 decorin -- Integra, Telios
 defensins -- Large Scale Biology
 DEGR-VIIa
 Delimmunised antibody 3B6/22 AGEN
 Deimmunised anti-cancer antibodies --
 Biovation/Viragen
 Dendroamide A
 Dengue vaccine -- Bavarian Nordic, Merck
 denileukin diftotox -- Ligand
 DES-1101 -- Desmos
 desirudin -- Novartis
 desmopressin -- Unigene
 Desmoteplase -- Merck, Schering AG
 Destablase
 Diabetes gene therapy -- DeveloGen, Pfizer
 Diabetes therapy -- Crucell
 Diabetes type 1 vaccine -- Diamyd
 Therapeutics
 DiaCIM -- YM BioSciences
 dialytic oligopeptides -- Research Corp
 Diamyd -- Diamyd Therapeutics
 DiaPep227 -- Peppen
 DiavaX -- Corixa
 Digoxin MAb -- Glaxo
 Diphtheria tetanus pertussis-hepatitis B
 vaccine -- GlaxoSmithKline
 DIR therapy -- Solis Therapeutics --
 DNase -- Genentech
 Dornase alfa -- Genentech
 Dornase alfa, inhalation -- Genentech
 Doxorubicin-anti-CEA MAb conjugate --
 Immunomedics
 DP-107 -- Trimeris
 drotrecogin alfa -- Eli Lilly
 DTctGMCSF
 DTP-polio vaccine -- Aventis Pasteur
 DU 257-KM231 antibody conjugate --
 Kyowa
 dural graft matrix -- Integra
 Duteplase -- Baxter Intl.
 DWP-401 -- Daewoong
 DWP-404 -- Daewoong
 DWP-408 -- Daewoong
 Dx 88 (Epi-KAL2) -- Dyax
 Dx 890 (elastin inhibitors) -- Dyax
 E coli O157 vaccine -- NIH
 E21-R -- BresaGen
 Eastern equine encephalitis virus vaccine --
 Echicetin --
 Echinhibin 1 --
 Echistatin -- Merck
 Echitamine --
 Ecromeximab -- Kyowa Hakko
 EC-SOD -- PPL Therapeutics
 Eculizumab (5G1.1) -- Alexion
 EDF -- Ajinomoto
 EDN derivative -- NIH
 EDNA -- NIH
 Edobacomab -- XOMA
 Edrecolomab -- Centocor
 EF 5077
 Efalizumab -- Genentech
 EGF fusion toxin -- Seragen, Ligand
 EGF-P64k vaccine -- Center of Molecular
 Immunology
 EL 246 -- LigoCyte
 elastase inhibitor -- Synergen
 elcatonin -- Therapicon
 EMD 72000 -- Merck KGaA
 Emdogain -- BIORA
 emfilermin -- AMRAD
 Emoctakin -- Novartis
 enamel matrix protein -- BIORA
 Endo III -- NYU
 endostatin -- EntreMed, Pharis
 Enhancins -- Micrologix
 Enlimomab -- Isis Pharm.
 Enoxaparin sodium -- Phamuka

FIG. 28J

41/498

enzyme linked antibody nutrient depletion
 therapy -- KS Biomedix Holdings
 Eosinophil-derived neutralizing agent --
 EP-51216 -- Asta Medica
 EP-51389 -- Asta Medica
 EPH family ligands -- Regeneron
 Epidermal growth factor -- Hitachi Kasei,
 Johnson & Johnson
 Epidermal growth factor fusion toxin --
 Senetek
 Epidermal growth factor-genistein --
 EPI-HNE-4 -- Dyax
 EPI-KAL2 -- Dyax
 Epoetin-alfa -- Amgen, Dragon
 Pharmaceuticals, Nanjing Huaxin
 Epratuzumab -- Immunomedics
 Epstein-Barr virus vaccine --
 Aviron/SmithKline Beecham, Bioresearch
 Eptacog alfa -- Novo Nordisk
 Eptifibatide -- COR Therapeutics
 erb-38 --
 Erlizumab -- Genentech
 erythropoietin -- Alkermes, ProLease, Dong-
 A, Elanex, Genetics Institute, LG Chem,
 Protein Sciences, Serono, Snow Brand,
 SRC VB VECTOR, Transkaryotic
 Therapies
 Erythropoietin Beta -- Hoffman La Roche
 Erythropoietin/Epoetin alfa -- Chugai
 Escherichia coli vaccine -- North American
 Vaccine, SBL Vaccin, Swiss Serum and
 Vaccine Institute Berne
 etanercept -- Immunex
 examorelin -- Mediolanum
 Exendin 4 -- Amylin
 exonuclease VII
 F 105 -- Centocor
 F-992 -- Fornix
 Factor IX -- Alpha Therapeutics, Welfide
 Corp., CSL, enetics Institute/AHP,
 Pharmacia, PPL Therapeutics
 Factor IX gene therapy -- Cell Genesys
 Factor VII -- Novo Nordisk, Bayer, Baxter
 Intl.
 Factor VIIa -- PPL Therapeutics,
 ZymoGenetics
 Factor VIII -- Bayer Genentech, Beaufour-
 Ipsen, CLB, Inex, Octagen, Pharmacia,
 Pharming
 Factor VIII -- PEGylated -- Bayer
 Factor VIII fragments -- Pharmacia
 Factor VIII gene therapy -- Targeted
 Genetics
 Factor VIII sucrose formulation -- Bayer,
 Genentech
 Factor VIII-2 -- Bayer
 Factor VIII-3 -- Bayer
 Factor Xa inhibitors -- Merck, Novo Nordisk,
 Mochida
 Factor XIII -- ZymoGenetics
 Factors VIII and IX gene therapy -- Genetics
 Institute/Targeted Genetics
 Famoxin -- Genset
 Fas (delta) TM protein -- LXR BioTech.
 Fas TR -- Human Genome Sciences
 Felvizumab -- Scotgen
 FFR-VIIa -- Novo Nordisk
 FG-001 -- F-Gene
 FG-002 -- F-Gene
 FG-004 -- F-Gene
 FG-005 -- F-Gene
 FGF + fibrin -- Repair
 Fibrimage -- Bio-Tech. General
 fibrin-binding peptides -- ISIS Innovation
 fibrinogen -- PPL Therapeutics, Pharming
 fibroblast growth factor -- Chiron, NYU,
 Ramot, ZymoGenetics
 fibrolase conjugate -- Schering AG
 Filgrastim -- Amgen
 filgrastim -- PDA modified -- Xencor
 FLT-3 ligand -- Immunex
 FN18 CRM9 --

FIG. 28K

42/498

follistatin -- Biotech Australia, Human Therapeutics
 follitropin alfa -- Alkermes, ProLease, PowderJect, Serono, Akzo Nobel
 Follitropin Beta -- Bayer, Organon
 FP 59
 FSH -- Ferring
 FSH + LH -- Ferring
 F-spondin -- CeNeS
 fusion protein delivery system -- UAB Research Foundation
 fusion toxins -- Boston Life Sciences
 G 5598 -- Genentech
 GA-II -- Transkaryotic Therapies
 Gamma-interferon analogues -- SRC VB VECTOR
 Ganirelix -- Roche
 gastric lipase -- Meristem
 Gavilimomab --
 G-CSF -- Amgen, SRC VB VECTOR
 GDF-1 -- CeNeS
 GDF-5 -- Biopharm
 GDNF (glial derived neurotrophic factor) -- Amgen
 gelsolin -- Biogen
 Gemtuzumab ozogamicin -- Celltech
 Gene-activated epoetin-alfa -- Aventis Pharma -- Transkaryotic Therapies
 Glanzmann thrombasthenia gene therapy --
 Glatiramer acetate -- Yeda
 glial growth factor 2 -- CeNeS
 GLP-1 -- Amylin, Suntory, TheraTech, Watson
 GLP-1 peptide analogues -- Zealand Pharmaceuticals
 GLP-2 -- Novo Nordisk, Ontario, Inc., Suntory Limited
 glucagon -- Eli Lilly, ZymoGenetics
 Glucagon-like peptide 1-7-36 amide -- Suntory
 Glucogen-like peptide -- Amylin
 Glucocerebrosidase -- Genzyme
 glutamate decarboxylase -- Genzyme Transgenics
 Glycoprotein S3 -- Kureha
 GM-CSF -- Immuhex
 GM-CSF tumour vaccine -- PowderJect
 GnRH immunotherapeutic -- Protherics
 Goserelin (LhRH antagonist) -- AstraZeneca
 gp75 antigen -- ImClone
 gp96 -- Antigenics
 GPI 0100 -- Galenica
 GR 4991W93 -- GlaxoSmithKline
 Granulocyte colony-stimulating factor -- Dong-A
 Granulocyte colony-stimulating factor conjugate
 grass allergy therapy -- Dynavax
 GRF1-44 -- ICN
 Growth Factor -- Chiron, Atrigel, Atrix, Innogenetics, ZymoGenetics, Novo
 growth factor peptides -- Biotherapeutics
 growth hormone -- LG Chem
 growth hormone, Recombinant human -- Serono
 GT 4086 -- Gliatech
 GW 353430 -- GlaxoSmithKline
 GW-278884 -- GlaxoSmithKline
 H 11 -- Viventia Biotech
 H5N1 influenza A virus vaccine -- Protein Sciences
 haemoglobin -- Biopure
 haemoglobin 3011, Recombinant -- Baxter Healthcare
 haemoglobin crosumaril -- Baxter Intl.
 haemoglobin stabilized -- Ajinomoto
 haemoglobin, recombinant -- Apex
 HAF -- Immune Response
 Hantavirus vaccine
 HB 19
 HBNF -- Regeneron
 HCC-1 -- Pharis

FIG. 28L

43/498

hCG -- Milkhaus
 hCG vaccine -- Zonagen
 HE-317 -- Hollis-Eden Pharmaceuticals
 Heat shock protein cancer and influenza
 vaccines -- StressGen
 Helicobacter pylori vaccine -- Acambis,
 AstraZeneca/CSL, Chiron, Provalis
 Helistat-G -- GalaGen
 Hemolink -- Hemosol
 hepapoietin -- Snow Brand
 heparanase -- InSight
 heparinase I -- Ibex
 heparinase III -- Ibex
 Hepatitis A vaccine -- American Biogenetic
 Sciences
 Hepatitis A vaccine inactivated
 Hepatitis A vaccine Nothav -- Chiron
 Hepatitis A-hepatitis B vaccine --
 GlaxoSmithKline
 hepatitis B therapy -- Tripep
 Hepatitis B vaccine -- Amgen, Chiron SpA,
 Meiji Milk, NIS, Prodeva, PowderJect,
 Rhein Biotech
 Hepatitis B vaccine recombinant -- Evans
 Vaccines, Epitex Combiotech, Genentech,
 MedImmune, Merck Sharp & Dohme,
 Rhein Biotech, Shantha Biotechnics,
 Vector, Yeda
 Hepatitis B vaccine recombinant TGP 943 --
 Takeda
 Hepatitis C vaccine -- Bavarian Nordic,
 Chiron, Innogenetics Acambis,
 Hepatitis D vaccine -- Chiron Vaccines
 Hepatitis E vaccine recombinant --
 Genelabs/GlaxoSmithKline, Novavax
 hepatocyte growth factor -- Panorama,
 Sosei
 hepatocyte growth factor kringle fragments -
 - EntreMed
 Her-2/Neu peptides -- Corixa
 Herpes simplex glycoprotein DNA vaccine --
 Merck, Wyeth-Lederle Vaccines-Malvern,
 Genentech, GlaxoSmithKline, Chiron,
 Takeda
 Herpes simplex vaccine -- Cantab
 Pharmaceuticals, CEL-SCI, Henderson
 Morley
 Herpes simplex vaccine live -- ImClone
 Systems/Wyeth-Lederle, Aventis Pasteur
 HGF derivatives -- Dompe
 hAPP vaccine -- Crucell
 Hib-hepatitis B vaccine -- Aventis Pasteur
 HIC 1
 HIP-- Altachem
 Hirudins -- Biopharma, Cangene, Dongkook,
 Japan Energy Corporation, Pharmacia
 Corporation, SIR International, Sanofi-
 Synthelabo, Sotragene, Rhein Biotech
 HIV edible vaccine -- ProdiGene
 HIV gp120 vaccine -- Chiron, Ajinomoto,
 GlaxoSmithKline, ID Vaccine, Progenics,
 VaxGen
 HIV gp120 vaccine gene therapy --
 HIV gp160 DNA vaccine -- PowderJect,
 Aventis Pasteur, Oncogen, Hyland
 Immuno, Protein Sciences
 HIV gp41 vaccine -- Panacos
 HIV HGP-30W vaccine -- CEL-SCI
 HIV immune globulin -- Abbott, Chiron
 HIV peptides -- American Home Products
 HIV vaccine -- Applied bioTech., Axis
 Genetics, Biogen, Bristol-Myers Squibb,
 Genentech, Korea Green Cross, NIS,
 Oncogen, Protein Sciences Corporation,
 Terumo, Tonen Corporation, Wyeth-
 Ayerst, Wyeth-Lederle Vaccines-Malvern,
 Advanced BioScience Laboratories,
 Bavarian Nordic, Bavarian Nordic/Statens
 Serum Institute, GeneCure, Immune
 Response, Progenics, Theron Biologics,
 United Biomedical, Chiron

FIG. 28M

44/498

HIV vaccine vCP1433 -- Aventis Pasteur
 HIV vaccine vCP1452 -- Aventis Pasteur
 HIV vaccine vCP205 -- Aventis Pasteur
 HL-9 -- American BioScience
 HM-9239 -- Cytran
 HML-103 -- Hemosol
 HML-104 -- Hemosol
 HML-105 -- Hemosol
 HML-109 -- Hemosol
 HML-110 -- Hemosol
 HML-121 -- Hemosol
 hNLP -- Pharis
 Hookworm vaccine
 host-vector vaccines -- Henogen
 HPM 1 -- Chugai
 HPV vaccine -- MediGene
 HSA -- Meristem
 HSF -- StressGen
 HSP carriers -- Weizmann, Yeda, Peptor
 HSPPC-70 -- Antigenics
 HSPPC-96, pathogen-derived -- Antigenics
 HSV 863 -- Novartis
 HTLV-I DNA vaccine
 HTLV-I vaccine
 HTLV-II vaccine -- Access
 HU 901 -- Tanox
 Hu23F2G -- ICOS
 HuHMF1
 HumaLYM -- Intracell
 Human krebs statika -- Yamanouchi
 human monoclonal antibodies --
 Abgenix/Biogen, Abgenix/ Corixa,
 Abgenix/immunex, Abgenix/Lexicon,
 Abgenix/ Pfizer, Athersys/Medarex,
 Biogen/MorphoSys, CAT/Searle,
 Centocor/Medarex, Corixa/Kirin Brewery,
 Corixa/Medarex, Eos BioTech./Medarex,
 Eos/Xenex, Exelixis/Protein Design
 Labs, ImmunoGen/ Raven, Medarex/
 B.Twelve, MorphoSys/ImmunoGen, XTL
 Biopharmaceuticals/Dyax,
 Human monoclonal antibodies --
 Medarex/Northwest Biotherapeutics,
 Medarex/Seattle Genetics
 human netrin-1 -- Exelixis
 human papillomavirus antibodies -- Epicyte
 Human papillomavirus vaccine -- Biotech
 Australia, IDEC, StressGen
 Human papillomavirus vaccine MEDI 501 --
 MedImmune/GlaxoSmithKline
 Human papillomavirus vaccine MEDI
 503/MEDI 504 --
 MedImmune/GlaxoSmithKline
 Human papillomavirus vaccine TA-CIN --
 Cantab Pharmaceuticals
 Human papillomavirus vaccine TA-HPV --
 Cantab Pharmaceuticals
 Human papillomavirus vaccine TH-GW --
 Cantab/GlaxoSmithKline
 human polyclonal antibodies -- Biosite/Eos
 BioTech./ Medarex
 human type II anti factor VIII monoclonal
 antibodies -- ThromboGenics
 humanised anti glycoprotein Ib murine
 monoclonal antibodies -- ThromboGenics
 HumaRAD -- Intracell
 HuMax EGFR -- Genmab
 HuMax-CD4 -- Medarex
 HuMax-IL15 -- Genmab
 HYB 190 -- Hybridon
 HYB 676 -- Hybridon
 I-125 Mab A33 -- Celltech
 Ibritumomab tiuxetan -- IDEC
 IBT-9401 -- Ibx
 IBT-9402 -- Ibx
 IC 14 -- ICOS
 Idarubicin anti-Ly-2.1 --
 IDEC 114 -- IDEC
 IDEC 131 -- IDEC
 IDEC 152 -- IDEC
 IDM 1 -- IDM
 IDPS -- Hollis-Eden Pharmaceuticals

FIG. 28N

45/498

iduronate-2-sulfatase -- Transkaryotic Therapies
 IGF/IBP-2-13 -- Pharis
 IGN-101 -- Igeneon
 IK HIR02 -- Iketon
 IL-11 -- Genetics Institute/AHP
 IL-13-PE38 -- NeoPharm
 IL-17 receptor -- Immunex
 IL-18BP -- Yeda
 IL-1Hy1 -- Hyseq
 IL-1 β -- Celltech
 IL-1 β adjuvant -- Celltech
 IL-2 -- Chiron
 IL-2 + IL-12 -- Hoffman La-Roche
 IL-6/sIL-6R fusion -- Hadasit
 IL-6R derivative -- Tosoh
 IL-7-Dap 389 fusion toxin -- Ligand
 IL-21 -- Novo Nordisk, ZymoGenetics
 IM-862 -- Cytran
 IMC-1C11 -- ImClone
 imiglucerase -- Genzyme
 Immune globulin intravenous (human) -- Hoffman La Roche
 immune privilege factor -- Proneuron
 Immunocal -- Immunotec
 Immunogene therapy -- Briana Bio-Tech
 Immunoliposomal 5-fluorodeoxyuridine-dipalmitate --
 immunosuppressant vaccine -- Aixlie
 immunotoxin -- Antisoma, NIH
 ImmuRAIT-Re-188 -- Immunomedics
 imreg-1 -- Imreg
 infertility -- Johnson & Johnson, E-TRANS
 Infliximab -- Centocor
 Influenza virus vaccine -- Aventis Pasteur, Protein Sciences
 inhibin -- Biotech Australia, Human Therapeutics
 Inhibitory G protein gene therapy
 INKP-2001 -- InKine
 Inolimomab -- Diaclone
 insulin -- AutoImmune, Altea, Biobras, BioSante, Bio-Tech. General, Chong Kun Dang, Emisphere, Flamel, Provalis, Rhein Biotech, TranXenoGen
 insulin (bovine) -- Novartis
 insulin analogue -- Eli Lilly
 Insulin Aspart -- Novo Nordisk
 insulin detemir -- Novo Nordisk
 insulin glargine -- Aventis
 insulin inhaled -- Inhale Therapeutics Systems, Alkermes
 insulin oral -- Inovax
 insulin, AeroDose -- AeroGen
 insulin, AERx -- Aradigm
 insulin, BEODAS -- Elan
 insulin, Biphasix -- Helix
 insulin, buccal -- Generex
 insulin, I2R -- Flemington
 insulin, intranasal -- Bentley
 insulin, oral -- Nobex, Unigene
 insulin, Orasome -- Endorex
 insulin, ProMaxx -- Epic
 insulin, Quadrant -- Elan
 insulin, recombinant -- Aventis
 insulin, Spiros -- Elan
 insulin, Transfersome -- IDEA
 insulin, Zymo, recombinant -- Novo Nordisk
 insulinotropin -- Scios
 Insulysin gene therapy --
 integrin antagonists -- Merck
 interferon (Alpha2) -- SRC VB VECTOR, Viragen, Dong-A, Hoffman La-Roche, Genentech
 interferon -- BioMedicines, Human Genome Sciences
 interferon (Alfa-n3) -- Interferon Sciences Intl.
 interferon (Alpha), Biphasix -- Helix

FIG. 280

46/498

interferon (Alpha)—Amgen, BioNative,
 Novartis, Genzyme Transgenics,
 Hayashibara, Inhale Therapeutics
 Systems, Medusa, Flamel, Dong-A,
 GeneTrol, Nasteck, Shantha,
 Wassermann, LG Chem, Sumitomo,
 Aventis, Behring EGIS, Pepgen, Servier,
 Rhein Biotech,
 interferon (Alpha2A)
 interferon (Alpha2B) — Enzon, Schering-
 Plough, Biogen, IDEA
 interferon (Alpha-N1) — GlaxoSmithKline
 interferon (beta) — Rentschler, GeneTrol,
 Meristem, Rhein Biotech, Toray, Yeda,
 Daiichi, Mochida
 interferon (Beta1A) — Sero, Biogen
 interferon (beta1A), inhale — Biogen
 interferon (β1b)— Chiron
 interferon (tau)— Pepgen
 Interferon alfacon-1 — Amgen
 Interferon alpha-2a vaccine
 Interferon Beta 1b — Schering/Chiron,
 InterMune
 Interferon Gamma — Boehringer Ingelheim,
 Sheffield, Rentschler, Hayashibara
 interferon receptor, Type I — Sero
 interferon (Gamma1B) — Genentech
 Interferon-alpha-2b + ribavirin — Biogen,
 ICN
 Interferon-alpha-2b gene therapy —
 Schering-Plough
 Interferon-con1 gene therapy —
 interleukin-1 antagonists — Dompe
 Interleukin-1 receptor antagonist — Abbott
 Bioresearch, Pharmacia
 Interleukin-1 receptor type I — Immunex
 interleukin-1 receptor Type II — Immunex
 Interleukin-1 trap — Regeneron
 Interleukin-1-alpha — Immunex/Roche
 interleukin-2 — SRC VB VECTOR,
 Ajinomoto, Biomira, Chiron
 IL-2/ diphtheria toxin — Ligand
 Interleukin-3 — Cangene
 Interleukin-4 — Immunology Ventures,
 Sanofi Winthrop, Schering-Plough,
 Immunex/ Sanofi Winthrop, Bayer, Ono
 interleukin-4 + TNF-Alpha — NIH
 interleukin-4 agonist — Bayer
 interleukin-4 fusion toxin — Ligand
 Interleukin-4 receptor — Immunex, Immun
 Interleukin-6 — Ajinomoto, Cangene, Yeda,
 Genetics Institute, Novartis
 interleukin-6 fusion protein
 interleukin-6 fusion toxin — Ligand, Sero
 interleukin-7 — IC Innovations
 interleukin-7 receptor — Immunex
 interleukin-8 antagonists — Kyowa
 Hakko/Millennium/Pfizer
 interleukin-9 antagonists — Genaera
 Interleukin-10 — DNAX, Schering-Plough
 Interleukin-10 gene therapy —
 interleukin-12 — Genetics Institute, Hoffman
 La-Roche
 interleukin-13 — Sanofi
 interleukin-13 antagonists — AMRAD
 Interleukin-13-PE38QQR
 interleukin-15 — Immunex
 interleukin-16 — Research Corp
 Interleukin-18 — GlaxoSmithKline
 Interleukin-18 binding protein — Sero
 Ior-P3 — Center of Molecular Immunology
 IP-10 — NIH
 IPF — Metabolex
 IR-501 — Immune Response
 ISIS 9125 — Isis Pharmaceuticals
 ISURF No. 1554 — Millennium
 ISURF No. 1866 — Iowa State Univer.
 ITF-1697 — Italfarmaco
 IxC 162 — Ixion
 J 695 — Cambridge Antibody Tech.,
 Genetics Inst., Knoll
 Jagged + FGF — Repair

FIG. 28P

47/498

JKC-362 -- Phoenix Pharmaceuticals
 JTP-2942 -- Japan Tobacco
 Juman monoclonal antibodies --
 Medarex/Raven
 K02 -- Axyx Pharmaceuticals
 Keliximab -- IDEC
 Keyhole limpet haemocyanin
 KGF -- Amgen
 KM 871 -- Kyowa
 KPI 135 -- Scios
 KPI-022 -- Scios
 Kringle 5
 KSB 304
 KSB-201 -- KS Biomedex
 L 696418 -- Merck
 L 703801 -- Merck
 L1 -- Acorda
 L-761191 -- Merck
 lactoferrin -- Meristem, Pharming, Agennix
 lactoferrin cardio -- Pharming
 LAG-3 -- Seroxo
 LAIT -- GEMMA
 LAK cell cytotoxin -- Arizona
 lamellarins -- PharmaMar/University of
 Malaga
 laminin A peptides -- NIH
 lanoteplase -- Genetics Institute
 laronidase -- BioMarin
 Lassa fever vaccine
 LCAT -- NIH
 LDP 01 -- Millennium
 LDP 02 -- Millennium
 Lecithinized superoxide dismutase --
 Seikagaku
 LeIF adjuvant -- Corixa
 leishmaniasis vaccine -- Corixa
 lenercept -- Hoffman La-Roche
 Lenograstim -- Aventis, Chugai
 lepirudin -- Aventis
 leptin -- Amgen, IC Innovations
 Leptin gene therapy -- Chiron Corporation
 leptin, 2nd-generation -- Amgen
 leridistim -- Pharmacia
 leuprolide, ProMaxx -- Epic
 leuprorelin, oral -- Unigene
 LeuTech -- Papatin
 LEX 032 -- SuperGen
 LiDEPT -- Novartis
 Lintuzumab (anti-CD33 MAb) -- Protein
 Design Labs
 lipase -- Altus Biologics
 lipid A vaccine -- EntreMed
 lipid-linked anchor Tech. -- ICRT, ID
 Biomedical
 liposome-CD4 Tech. -- Sheffield
 Listeria monocytogenes vaccine
 LMB 1
 LMB 7
 LMB 9 -- Battelle Memorial Institute, NIH
 LM-CD45 -- Cantab Pharmaceuticals
 lovastatin -- Merck
 LSA-3
 LT- β receptor -- Biogen
 lung cancer vaccine -- Corixa
 lusupultide -- Scios
 L-Vax -- AVAX
 LY 355455 -- Eli Lilly
 LY 366405 -- Eli Lilly
 LY-355101 -- Eli Lilly
 Lyme disease DNA vaccine -- Vical/Aventis
 Pasteur
 Lyme disease vaccine -- Aquila
 Biopharmaceuticals, Aventis, Pasteur,
 Symbicom, GlaxoSmithKline, Hyland
 Immuno, MedImmune
 Lymphocytic choriomeningitis virus vaccine
 lymphoma vaccine -- Biomira, Genitope
 LYP18
 lys plasminogen, recombinant
 Lysosomal storage disease gene therapy --
 Avigen
 lysostaphin -- Nutrition 21

FIG. 28Q

48/498

M 23 -- Gruenenthal
 M1 monoclonal antibodies -- Acorda
 Therapeutics
 MA 16N7C2 -- Corvas Intl.
 malaria vaccine -- GlaxoSmithKline,
 AdProTech, Antigenics, Apovia, Aventis
 Pasteur, Axis Genetics, Behringwerke,
 CDCP, Chiron Vaccines, Genzyme
 Transgenics, Hawaii, MedImmune, NIH,
 NYU, Oxxon, Roche/Saramane, Biotech
 Australia, Rx Tech
 Malaria vaccine CDC/NIIMALVAC-1
 malaria vaccine, multicomponent
 mammaglobin -- Corixa
 mammastatin -- Biotherapeutics
 mannan-binding lectin -- NatlImmu
 mannan-MUC1 -- Psiron
 MAP 30
 Marinovir -- Phytera
 MARstem -- Maret
 MB-015 -- Mochida
 MBP -- ImmuLogic
 MCI-028 -- Mitsubishi-Tokyo
 MCIF -- Human Genome Sciences
 MDC -- Advanced BioScience -- Akzo
 Nobel, ICOS
 MDX 11 -- Medarex
 MDX 210 -- Medarex
 MDX 22 -- Medarex
 MDX 22
 MDX 240 -- Medarex
 MDX 33
 MDX 44 -- Medarex
 MDX 447 -- Medarex
 MDX H210 -- Medarex
 MDX RA -- Houston BioTech., Medarex
 ME-104 -- Pharmexa
 Measles vaccine
 Mecasermin -- Cephalon/Chiron, Chiron
 MEDI 488 -- MedImmune
 MEDI 500
 MEDI 507 -- BioTransplant
 melanin concentrating hormone --
 Neurocrine Biosciences
 melanocortins -- OMRF
 Melanoma monoclonal antibodies -- Viragen
 melanoma vaccine -- GlaxoSmithKline,
 Akzo Nobel, Avant, Aventis Pasteur,
 Bavarian Nordic, Biovector, CancerVax,
 Genzyme Molecular Oncology, Humbolt,
 ImClone Systems, Memorial, NYU, Oxxon
 Melanoma vaccine Magevac -- Therion
 memory enhancers -- Scios
 meningococcal B vaccine -- Chiron
 meningococcal vaccine -- CAMR
 Meningococcal vaccine group B conjugate -
 - North American Vaccine
 Meningococcal vaccine group B
 recombinant -- BioChem Vaccines,
 Microscience
 Meningococcal vaccine group Y conjugate -
 - North American Vaccine
 Meningococcal vaccine groups A B and C
 conjugate -- North American Vaccine
 Mepolizumab -- GlaxoSmithKline
 Metastatin -- EntreMed, Takeda
 Met-CkB7 -- Human Genome Sciences
 met-enkephalin -- TNI
 METH-1 -- Human Genome Sciences
 methioninase -- AntiCancer
 Methionine lyase gene therapy --
 AntiCancer
 Met-RANTES -- Genexa Biomedical,
 Seroxo
 Metreleptin
 Microtubule inhibitor MAb
 Immunogen/Abgenix
 MGDF -- Kirin
 MGv -- Progenics
 micrin -- Endocrine
 microplasmin -- ThromboGenics
 MIF -- Genetics Institute

FIG. 28R

49/498

migration inhibitory factor -- NIH	MAb 45-2D9- -- haematoporphyrin conjugate
Mim CD4.1 -- Xycte Therapies	MAb 4B4
mirostipen -- Human Genome Sciences	MAb 4E3-CPA conjugate -- BCM Oncologia
Mitumomab (BEC-2) -- ImClone Systems, Merck KGaA	MAb 4E3-daunorubicin conjugate
MK 852 -- Merck	MAb 50-6
MLN 1202 (Anti-CCR2 monoclonal antibody) -- Millenium Pharmaceuticals	MAb 50-61A -- Institut Pasteur
Mobenakin -- NIS	MAb 5A8 -- Biogen
molgramostim -- Genetics Institute, Novartis	MAb 791T/36-methotrexate conjugate
monoclonal antibodies -- Abgenix/Celltech, Immusol/ Medarex, Viragen/ Roslin Institute, Cambridge Antibody Tech./Elan	MAb 7c11.e8
MAb 108 --	MAb 7E11 C5-selenocystamine conjugate
MAb 10D5 --	MAb 93KA9 -- Novartis
MAb 14.18-interleukin-2 immunocytokine -- Lexigen	MAb A5B7-cisplatin conjugate -- Biodynamics Research, Pharmacia
MAb 14G2a --	MAb A5B7-I-131
MAb 15A10 --	MAb A7
MAb 170 -- Biomira	MAb A717 -- Exocell
MAb 177Lu CC49 --	MAb A7-zinostatin conjugate
MAb 17F9	MAb ABX-RB2 -- Abgenix
MAb 1D7	MAb ACA 11
MAb 1F7 -- Immune Network	MAb AFP-I-131 -- Immunomedics
MAb 1H10-doxorubicin conjugate	MAb AP1
MAb 26-2F	MAb AZ1
MAb 2A11	MAb B3-LysPE40 conjugate
MAb 2E1 -- RW Johnson	MAb B4 -- United Biomedical
MAb 2F5	MAb B43 Genistein-conjugate
MAb 31.1 -- International BioImmune Systems	MAb B43.13-Tc-99m -- Biomira
MAb 32 -- Cambridge Antibody Tech., Peptech	MAb B43-PAP conjugate
MAb 323A3 -- Centocor	MAb B4G7-gelonin conjugate
MAb 3C5	MAb BCM 43-daunorubicin conjugate -- BCM Oncologia
MAb 3F12	MAb BIS-1
MAb 3F8	MAb BMS 181170 -- Bristol-Myers Squibb
MAb 42/6	MAb BR55-2
MAb 425 -- Merck KGaA	MAb BW494
MAb 447-52D -- Merck Sharp & Dohme	MAb C 242-DM1 conjugate -- ImmunoGen
	MAb C242-PE conjugate
	MAb c30-6
	MAb CA208-cytorhodin-S conjugate -- Hoechst Japan
	MAb CC49 -- Enzo

50/498

MAb ch14.18 --	MAb LL2-I-131 -- Immunomedics
MAb CH14.18-GM-CSF fusion protein --	MAb LL2-Y-90
Lexigen	MAb LS2D617 -- Hybritech
MAb chCE7	MAb LYM-1-gelonin conjugate
MAb CI-137 -- AMRAD	MAb LYM-1-I-131
MAb cisplatin conjugate	MAb LYM-1-Y-90
MAb CLB-CD19	MAb LYM-2 -- Peregrine
MAb CLB-CD19v	MAb M195
MAb CLL-1 -- Peregrine	MAb M195-bismuth 213 conjugate --
MAb CLL-1-GM-CSF conjugate	Protein Design Labs
MAb CLL-1-IL-2 conjugate -- Peregrine	MAb M195-gelonin conjugate
MAb CLN IgG -- doxorubicin conjugates	MAb M195-I-131
MAb conjugates -- Tanox	MAb M195-Y-90
MAb D612	MAb MA 33H1 -- Sanofi
MAb Dal B02	MAb MAD11
MAb DC101 -- ImClone	MAb MGb2
MAb EA 1 --	MAb MINT5
MAb EC708 -- Biovation	MAb MK2-23
MAb EP-5C7 -- Protein Design Labs	MAb MOC31 ETA(252-613) conjugate
MAb ERIC-1 -- ICRT	MAb MOC-31-In-111
MAb F105 gene therapy	MAb MOC-31-PE conjugate
MAb FC 2.15	MAb MR6 --
MAb G250 -- Centocor	MAb MRK-16 -- Aventis Pasteur
MAb GA6	MAb MS11G6
MAb GA733	MAb MX-DTPA BrE-3
MAb Gliomab-H -- Viventia Biotech	MAb MY9
MAb HB2-saporin conjugate	MAb Nd2 -- Tosoh
MAb HD 37 --	MAb NG-1 -- Hygeia
MAb HD37-ricin chain-A conjugate	MAb NM01 -- Nissin Food
MAb HNK20 -- Acambis	MAb OC 125
MAb huN901-DM1 conjugate --	MAb OC 125-CMA conjugate
ImmunoGen	MAb OKI-1 -- Ortho-McNeil
MAb I-131 CC49 -- Corixa	MAb OX52 -- Bioproducts for Science
MAb ICO25	MAb PMA5
MAb ICR12-CPG2 conjugate	MAb PR1
MAb ICR-62	MAb prost 30
MAb IRac-ricin A conjugate	MAb R-24
MAb K1	MAb R-24 α Human GD3 -- Celltech
MAb KS1-4-methotrexate conjugate	MAb RFB4-ricin chain A conjugate
MAb L6 -- Bristol-Myers Squibb, Oncogen	MAb RFT5-ricin chain A conjugate
MAb LICO 16-88	MAb SC 1

FIG. 28T

51/498

MAb SM-3 -- ICRT
 MAb SMART 1D10 -- Protein Design Labs
 MAb SMART ABL 364 -- Novartis
 MAb SN6f
 MAb SN6f-deglycosylated ricin A chain conjugate --
 MAb SN6j
 MAb SN7-ricin chain A conjugate
 MAb T101-Y-90 conjugate -- Hybritech
 MAb T-88 -- Chiron
 MAb TB94 -- Cancer Immunobiology
 MAb TEC 11
 MAb TES-23 -- Chugai
 MAb TM31 -- Avant
 MAb TNT-1 -- Cambridge Antibody Tech., Peregrine
 MAb TNT-3
 MAb TNT-3 -- IL2 fusion protein --
 MAb TP3-At-211
 MAb TP3-PAP conjugate --
 MAb UJ13A -- ICRT
 MAb UN3
 MAb ZME-018-gelonin conjugate
 MAb-BC2 -- GlaxoSmithKline
 MAb-DM1 conjugate -- ImmunoGen
 MAb-ricin-chain-A conjugate -- XOMA
 MAb-temoporfin conjugates
 Monopharm C -- Viventia Biotech
 montepelase -- Eisai
 montirelin hydrate -- Gruenenthal
 morotocog alfa -- Genetics Institute
 Morotocog-alfa -- Pharmacia
 MP 4
 MP-121 -- Biopharm
 MP-52 -- Biopharm
 MRA -- Chugai
 MS 28168 -- Mitsui Chemicals, Nihon Schering
 MSH fusion toxin -- Ligand
 MSI-99 -- Genæra
 MT 201 -- Micromet
 Muc-1 vaccine -- Corixa
 mucosal tolerance -- Aberdeen
 mullerian inhibiting subst
 muplestim -- Genetics Institute, Novartis,
 DSM Anti-Infectives
 murine MAb -- KS Biomedix
 Mutant somatropin -- JCR Pharmaceutical
 MV 833 -- Toagosei
 Mycoplasma pulmonis vaccine
 Mycoprex -- XOMA
 myeloperoxidase -- Henogen
 myostatin -- Genetics Institute
 Nacolomab tafenatox -- Pharmacia
 Nagrecor -- Scios
 nagrestipen -- British Biotech
 NAP-5 -- Corvas Intl.
 NAPc2 -- Corvas Intl.
 nartograstim -- Kyowa
 Natalizumab -- Protein Design Labs
 Nateplase -- NIH, Nihon Schering
 nateplase -- Schering AG
 NBI-3001 -- Neurocrine Biosci.
 NBI-5788 -- Neurocrine Biosci.
 NBI-6024 -- Neurocrine Biosci.
 Nef inhibitors -- BRI
 Neisseria gonorrhoea vaccine -- Antex Biologics
 Neomycin B-arginine conjugate
 Nerelimomab -- Chiron
 Nerve growth factor -- Amgen -- Chiron, Genentech
 Nerve growth factor gene therapy
 nesiritide citrate -- Scios
 neuregulin-2 -- CeNeS
 neurocan -- NYU
 neuronal delivery system -- CAMR
 Neurophil inhibitory Factor -- Corvas
 Neuroprotective vaccine -- University of Auckland
 neurotrophic chimaeras -- Regeneron
 neurotrophic factor -- NsGene, CereMedix

FIG. 28U

52/498

NeuroVax -- Immune Response
 neurturnin -- Genentech
 neutral endopeptidase -- Genentech
 NGF enhancers -- NeuroSearch
 NHL vaccine -- Large Scale Biology
 NIP45 -- Boston Life Sciences
 NKI-B20
 NM 01 -- Nissin Food
 NMI-139 -- NitroMed
 NMMP -- Genetics Institute
 NN-2211 -- Novo Nordisk
 Noggin -- Regeneron
 Nonacog alfa
 Norelin -- Biostar
 Norwalk virus vaccine
 NRLU 10 -- NeoRx
 NRLU 10 PE -- NeoRx
 NT-3 -- Regeneron
 NT-4/5 -- Genentech
 NU 3056
 NU 3076
 NX 1838 -- Gilead Sciences
 NY ESO-1/CAG-3 antigen -- NIH
 NYVAC-7 -- Aventis Pasteur
 NZ-1002 -- Novazyme
 obesity therapy -- Nobex
 OC 10426 -- Ontogen
 OC 144093 -- Ontogen
 OCIF -- Sankyo
 Oct-43 -- Otsuka
 Odulimomab -- Immunotech
 OK PSA - liposomal
 OKT3-gamma-1-ala-ala
 OM 991
 OM 992
 Omalizumab -- Genentech
 oncoimmunin-L -- NIH
 Oncolysin B -- ImmunoGen
 Oncolysin CD6 -- ImmunoGen
 Oncolysin M -- ImmunoGen
 Oncolysin S -- ImmunoGen
 Oncophage -- Antigenics
 Oncostatin M -- Bristol-Myers Squibb
 OncoVax-CL -- Jenner Biotherapies
 OncoVax-P -- Jenner Biotherapies
 onercept -- Yeda
 onychomycosis vaccine -- Boehringer
 Ingelheim
 opebecan -- XOMA
 opioids -- Arizona
 Oprelvekin -- Genetics Institute
 Oregovomab -- AltaRex
 Org-33408 b-- Akzo Nobel
 Orlip DP -- EpiCept
 oryzacystatin
 OSA peptides -- GenSci Regeneration
 osteoblast-cadherin GF -- Pharis
 Osteocalcin-thymidine kinase gene therapy
 osteogenic protein -- Curis
 osteopontin -- OraPharma
 osteoporosis peptides -- Integra, Telios
 osteoprotegerin -- Amgen, SnowBrand
 otitis media vaccines -- Antex Biologics
 ovarian cancer -- University of Alabama
 OX40-IgG fusion protein -- Cantab, Xenova
 P 246 -- Diatide
 P 30 -- Alfacell
 p1025 -- Active Biotech
 P-113^A -- Demegen
 P-16 peptide -- Transition Therapeutics
 p43 -- Ramot
 P-50 peptide -- Transition Therapeutics
 p53 + RAS vaccine -- NIH, NCI
 PACAP(1-27) analogue
 paediatric vaccines -- Chiron
 Pafase -- ICOS
 PAGE-4 plasmid DNA -- IDEC
 PAI-2 -- Biotech Australia, Human
 Therapeutics
 Palifermin (keratinocyte growth factor) --
 Amgen
 Palivizumab -- MedImmune

FIG. 28V

53/498

PAM 4 -- Merck
 pamiteplase -- Yamanouchi
 pancreatin, Minitabs -- Eurand
 Pangen -- Fournier
 Pantarin -- Selective Genetics
 Parainfluenza virus vaccine -- Pharmacia,
 Pierre Fabre
 paraoxanase -- Esperion
 parathyroid hormone -- Abiogen, Korea
 Green Cross
 Parathyroid hormone (1-34) --
 Chugai/Suntory
 Parkinson's disease gene therapy -- Cell
 Genesys/ Ceregene
 Parvovirus vaccine -- MedImmune
 PCP-Scan -- Immunomedics
 PDGF -- Chiron
 PDGF cocktail -- Theratechnologies
 peanut allergy therapy -- Dynavax
 PEG anti-ICAM MAb -- Boehringer
 Ingelheim
 PEG asparaginase -- Enzon
 PEG glucocerebrosidase
 PEG hirudin -- Knoll
 PEG interferon-alpha-2a -- Roche
 PEG interferon-alpha-2b + ribavirin --
 Biogen, Enzon, ICN Pharmaceuticals,
 Schering-Plough
 PEG MAb A5B7 --
 Pegacaristim -- Amgen -- Kirin Brewery --
 ZymoGenetics
 Pegaldesleukin -- Research Corp
 pegaspargase -- Enzon
 pegfilgrastim -- Amgen
 PEG-interferon Alpha -- Viragen
 PEG-interferon Alpha 2A -- Hoffman La-
 Roche
 PEG-interferon Alpha 2B -- Schering-
 Plough
 PEG-r-hirudin -- Abbott
 PEG-rHuMGDF -- Amgen
 PEG-uricase -- Mountain View
 Pegvisomant -- Genentech
 PEGylated proteins, PolyMASC -- Valentis
 PEGylated recombinant native human leptin
 -- Roche
 Pentumomab
 Penetratin -- Cyclacel
 Pepscan -- Antisoma
 peptide G -- Peptech, ICRT
 peptide vaccine -- NIH ,NCI
 Pexelizumab
 pexiganan acetate -- Genaera
 Pharmaprojects No. 3179 -- NYU
 Pharmaprojects No. 3390 -- Ernest Orlando
 Pharmaprojects No. 3417 -- Sumitomo
 Pharmaprojects No. 3777 -- Acambis
 Pharmaprojects No. 4209 -- XOMA
 Pharmaprojects No. 4349 -- Baxter Intl.
 Pharmaprojects No. 4651
 Pharmaprojects No. 4915 -- Avanir
 Pharmaprojects No. 5156 -- Rhizogenics
 Pharmaprojects No. 5200 -- Pfizer
 Pharmaprojects No. 5215 -- Origene
 Pharmaprojects No. 5216 -- Origene
 Pharmaprojects No. 5218 -- Origene
 Pharmaprojects No. 5267 -- ML
 Laboratories
 Pharmaprojects No. 5373 -- MorphoSys
 Pharmaprojects No. 5493 -- Metabolex
 Pharmaprojects No. 5707 -- Genentech
 Pharmaprojects No. 5728 -- Autogen
 Pharmaprojects No. 5733 -- BioMarin
 Pharmaprojects No. 5757 -- NIH
 Pharmaprojects No. 5765 -- Gryphon
 Pharmaprojects No. 5830 -- AntiCancer
 Pharmaprojects No. 5839 -- Dyax
 Pharmaprojects No. 5849 -- Johnson &
 Johnson
 Pharmaprojects No. 5860 -- Mitsubishi-
 Tokyo

FIG. 28W

54/498

- Pharmaprojects No. 5869 -- Oxford GlycoSciences
- Pharmaprojects No. 5883 -- Asahi Brewery
- Pharmaprojects No. 5947 -- StressGen
- Pharmaprojects No. 5961 -- Theratechnologies
- Pharmaprojects No. 5962 -- NIH
- Pharmaprojects No. 5966 -- NIH
- Pharmaprojects No. 5994 -- Pharming
- Pharmaprojects No. 5995 -- Pharming
- Pharmaprojects No. 6023 -- IMMUCON
- Pharmaprojects No. 6063 -- Cytoclonal
- Pharmaprojects No. 6073 -- SIDDCO
- Pharmaprojects No. 6115 -- Genzyme
- Pharmaprojects No. 6227 -- NIH
- Pharmaprojects No. 6230 -- NIH
- Pharmaprojects No. 6236 -- NIH
- Pharmaprojects No. 6243 -- NIH
- Pharmaprojects No. 6244 -- NIH
- Pharmaprojects No. 6281 -- Senetek
- Pharmaprojects No. 6365 -- NIH
- Pharmaprojects No. 6368 -- NIH
- Pharmaprojects No. 6373 -- NIH
- Pharmaprojects No. 6408 -- Pan Pacific
- Pharmaprojects No. 6410 -- Athersys
- Pharmaprojects No. 6421 -- Oxford GlycoSciences
- Pharmaprojects No. 6522 -- Maxygen
- Pharmaprojects No. 6523 -- Pharis
- Pharmaprojects No. 6538 -- Maxygen
- Pharmaprojects No. 6554 -- APALEXO
- Pharmaprojects No. 6560 -- Ardana
- Pharmaprojects No. 6562 -- Bayer
- Pharmaprojects No. 6569 -- Eos
- Phenoxazine
- Phenylase -- Ibbex
- Pigment epithelium derived factor -- plasminogen activator inhibitor-1, recombinant -- DuPont Pharmaceuticals
- Plasminogen activators -- Abbott Laboratories, American Home Products, Boehringer Mannheim, Chiron Corporation, DuPont Pharmaceuticals, Eli Lilly, Shionogi, Genentech, Genetics Institute, GlaxoSmithKline, Hemispherx Biopharma, Merck & Co, Novartis, Pharmacia Corporation, Wakamoto, Yeda
- plasminogen-related peptides -- Bio-Tech. General/MGH
- platelet factor 4 -- RepliGen
- Platelet-derived growth factor -- Amgen -- ZymoGenetics
- plusonemin -- Hayashibara
- PMD-2850 -- Protherics
- Pneumococcal vaccine -- Antex Biologics, Aventis Pasteur
- Pneumococcal vaccine intranasal -- BioChem Vaccines/Biovector
- PR1A3
- PR-39
- pralmorelin -- Kaken
- Pretarget-Lymphoma -- NeoRx
- Priliximab -- Centocor
- PRO 140 -- Progenics
- PRO 2000 -- Procept
- PRO 367 -- Progenics
- PRO 542 -- Progenics
- pro-Apo A-I -- Esperion
- prolactin -- Genzyme
- Prosaptide TX14(A) -- Bio-Tech. General
- prostate cancer antibodies -- Immunex, UroCor
- prostate cancer antibody therapy -- Genentech/UroGenesys, Genotherapeutics
- prostate cancer immunotherapeutics -- The PSMA Development Company
- prostate cancer vaccine -- Aventis Pasteur, Zonagen, Corixa, Dendreon, Jenner
- Biotherapies, Therion Biologics

FIG. 28X

55/498

prostate-specific antigen -- Entremed
 protein A -- RepliGen
 protein adhesives -- Enzon
 protein C -- Baxter Intl., PPL Therapeutics,
 ZymoGenetics
 protein C activator -- Gilead Sciences
 protein kinase R antagonists -- NIH
 protirelin -- Takeda
 protocadherin 2 -- Caprion
 Pro-urokinase -- Abbott, Bristol-Myers
 Squibb, Dainippon, Tosoh -- Welfide
 P-selectin glycoprotein ligand-1 -- Genetics
 Institute
 pseudomonal infections -- InterMune
 Pseudomonas vaccine -- CytoVax
 PSGL-Ig -- American Home Products
 PSP-94 -- Procyon
 PTH 1-34 -- Nobex
 Quilimmune-M -- Antigenics
 R 744 -- Roche
 R 101933
 R 125224 -- Sankyo
 RA therapy -- Cardion
 Rabies vaccine recombinant -- Aventis
 Pasteur, BioChem Vaccines, Kaketsuken
 Pharmaceuticals
 RadioTheraCIM -- YM BioSciences
 Ramot project No. 1315 -- Ramot
 Ramot project No. K-734A -- Ramot
 Ramot project No. K-734B -- Ramot
 Ranibizumab (Anti-VEGF fragment) --
 Genentech
 RANK -- Immunex
 ranpirinase -- Alfacell
 ranpirinase-anti-CD22 MAb -- Alfacell
 RANTES inhibitor -- Milan
 RAPID drug delivery systems -- ARIAD
 rasburicase -- Sanofi
 rBPI-21, topical -- XOMA
 RC 529 -- Corixa
 rCFTR -- Genzyme Transgenics

RD 62198
 rDnase -- Genentech
 RDP-58 -- SangStat
 RecepTox-Fce -- Keryx
 RecepTox-GnRH -- Keryx, MTR
 Technologies
 RecepTox-MBP -- Keryx, MTR
 Technologies
 recFSH -- Akzo Nobel, Organon
 REGA 3G12
 Regavirumab -- Teijin
 relaxin -- Connetics Corp
 Renal cancer vaccine -- MacroPharm
 repifermin -- Human Genome Sciences
 Respiratory syncytial virus PFP-2 vaccine --
 Wyeth-Lederle
 Respiratory syncytial virus vaccine --
 GlaxoSmithKline, Pharmacia, Pierre Fabre
 Respiratory syncytial virus vaccine
 inactivated
 Respiratory syncytial virus-parainfluenza
 virus vaccine -- Aventis Pasteur,
 Pharmacia
 Reteplase -- Boehringer Mannheim,
 Hoffman La-Roche
 Retropep -- Retroscreen
 RFB4 (dsFv) PE38
 RFI 641 -- American Home Products
 RFTS -- UAB Research Foundation
 RG 12986 -- Aventis Pasteur
 RG 83852 -- Aventis Pasteur
 RG-1059 -- RepliGen
 rGCR -- NIH
 rGLP-1 -- Restoragen
 rGRF -- Restoragen
 rh Insulin -- Eli Lilly
 RHAMM targeting peptides -- Cangene
 rHb1.1 -- Baxter Intl.
 rhCC10 -- Claragen
 rhCG -- SeroNo
 Rheumatoid arthritis gene therapy

FIG. 28Y

56/498

Rheumatoid arthritis vaccine -- Veterans

Affairs Medical Center

rhLH -- SeroNo

Ribozyme gene therapy -- Genset

Rickettsial vaccine recombinant

RIGScan CR -- Neoprobe

RIP-3 -- Rigel

Rituximab -- Genentech

RK-0202 -- RxKinetix

RLT peptide -- Esperion

rM/NEI -- IVAX

rmCRP -- Immtech

RN-1001 -- Renovo

RN-3 -- Renovo

RNase conjugate -- Immunomedics

RO 631908 -- Roche

Rotavirus vaccine -- Merck

RP 431 -- DuPont Pharmaceuticals

RP-128 -- Resolution

RPE65 gene therapy --

RPR 110173 -- Aventis Pasteur

RPR 115135 -- Aventis Pasteur

RPR 116258A -- Aventis Pasteur

rPSGL-Ig -- American Home Products

r-SPC surfactant -- Byk Gulden

RSV antibody -- Medimmune

Ruplizumab -- Biogen

rV-HER-2/neu -- Therion Biologics

SA 1042 -- Sankyo

sacrosidase -- Orphan Medical

Sant 7

Sargramostim -- Immunex

saruplase -- Gruenenthal

Satumomab -- Cytogen

SB 1 -- COR Therapeutics

SB 207448 -- GlaxoSmithKline

SB 208651 -- GlaxoSmithKline

SB 240683 -- GlaxoSmithKline

SB 249415 -- GlaxoSmithKline

SB 249417 -- GlaxoSmithKline

SB 6 -- COR Therapeutics

SB RA 31012 --

SC 56929 -- Pharmacia

SCA binding proteins -- Curis, Enzon

scFv(14E1)-ETA Berlex Laboratories,

Schering AG

ScFv(FRP5)-ETA --

ScFv6C6-PE40 --

SCH 55700 -- Celltech

Schistosomiasis vaccine -- Glaxo

Wellcome/Medeva, Brazil

SCPF -- Advanced Tissue Sciences

scuPA-suPAR complex -- Hadasit

SD-9427 -- Pharmacia

SDF-1 -- Ono

SDZ 215918 -- Novartis

SDZ 280125 -- Novartis

SDZ 89104 -- Novartis

SDZ ABL 364 -- Novartis

SDZ MMA 383 -- Novartis

Secretin -- Ferring, Repligen

serine protease inhbs -- Pharis

sermorelin acetate -- SeroNo

SERP-1 -- Viron

sertenef -- Dainippon

serum albumin, Recombinant human --

Aventis Behring

serum-derived factor -- Hadasit

Sevirumab -- Novartis

SGN 14 -- Seattle Genetics

SGN 15 -- Seattle Genetics

SGN 17/19 -- Seattle Genetics

SGN 30 -- Seattle Genetics

SGN-10 -- Seattle Genetics

SGN-11 -- Seattle Genetics

SH 306 -- DuPont Pharmaceuticals

Shanvac-B -- Shantha

Shigella flexneri vaccine -- Avant, Acambis,

Novavax

Shigella sonnei vaccine --

sICAM-1 -- Boehringer Ingelheim

Silteplase -- Genzyme

FIG. 28Z

57/498

SIV vaccine -- Endocon, Institut Pasteur
 SK 896 -- Sanwa Kagaku Kenkyusho
 SK-827 -- Sanwa Kagaku Kenkyusho
 Skeletex -- CellFactors
 SKF 106160 -- GlaxoSmithKline
 S-nitroso-AR545C --
 SNTP -- Active Biotech
 somatomedin-1 -- GroPep, Mitsubishi-
 Tokyo, NIH
 somatomedin-1 carrier protein -- Insmed
 somatostatin -- Ferring
 Somatotropin/
 Human Growth Hormone -- Bio-Tech.
 General, Eli Lilly
 somatropin -- Bio-Tech. General, Alkermes,
 ProLease, Aventis Behring, Biovector,
 Cangene, Dong-A, Eli Lilly, Emisphere,
 Enact, Genentech, Genzyme Transgenics,
 Grandis/InfiMed, CSL, InfiMed, MacroMed,
 Novartis, Novo Nordisk, Pharmacia
 Serono, TranXenoGen
 somatropin derivative -- Schering AG
 somatropin, AIR -- Eli Lilly
 Somatropin, inhaled -- Eli Lilly/Alkermes
 somatropin, Kabi -- Pharmacia
 somatropin, Orasome -- Novo Nordisk
 Sonermin -- Daiippon Pharmaceutical
 SP(V5.2)C -- Supertek
 SPf66
 sphingomyelinase -- Genzyme
 SR 29001 -- Sanofi
 SR 41476 -- Sanofi
 SR-29001 -- Sanofi
 SS1(dsFV)-PE38 -- NeoPharm
 β 2 microglobulin -- Avidex
 β 2-microglobulin fusion proteins -- NIH
 β -amyloid peptides -- CeNeS
 β -defensin -- Pharis
 Staphylococcus aureus infections --
 Inhibibex/ZLB
 Staphylococcus aureus vaccine conjugate --
 Nabi
 Staphylococcus therapy -- Tripep
 Staphylokinase -- Biovation, Prothera,
 Thrombogenetics
 Streptococcal A vaccine -- M6
 Pharmaceuticals, North American Vaccine
 Streptococcal B vaccine -- Microscience
 Streptococcal B vaccine recombinant --
 Biochem Vaccines
 Streptococcus pyogenes vaccine
 STRL-33 -- NIH
 Subalin -- SRC VB VECTOR
 SUIIS -- United Biomedical
 SUIIS-LHRH -- United Biomedical
 SUN-E3001 -- Suntory
 super high affinity monoclonal antibodies --
 YM BioSciences
 Superoxide dismutase -- Chiron, Enzon,
 Ube Industries, Bio-Tech, Yeda
 superoxide dismutase-2 -- OXIS
 suppressin -- UAB Research Foundation
 SY-161-P5 -- ThromboGenics
 SY-162 -- ThromboGenics
 Systemic lupus erythematosus vaccine --
 MedClone/VivoRx
 T cell receptor peptides -- Xoma
 T cell receptor peptide vaccine
 T4N5 liposomes -- AGI Dermatics
 TACI, soluble -- ZymoGenetics
 targeted apoptosis -- Antisoma
 tasonermin -- Boehringer Ingelheim
 TASP
 TASP-V
 Tat peptide analogues -- NIH
 TBP I -- Yeda
 TBP II
 TBV25H -- NIH
 Tc 99m ior cea1 -- Center of Molecular
 Immunology
 Tc 99m P 748 -- Diatide

FIG. 28AA

58/498

Tc 99m votumab -- Intracell
 Tc-99m rh-Annexin V -- Theseus Imaging
 teceleukin -- Biogen
 tenecteplase -- Genentech
 Teriparatide -- Armour Pharmaceuticals,
 Asahi Kasei, Eli Lilly
 terlipressin -- Ferring
 testisin -- AMRAD
 Tetra fibrin -- Roche
 TFPI -- EntrelMed
 tgD-IL-2 -- Takeda
 TGF-Alpha -- ZymoGenetics
 TGF- β -- Kolon
 TGF- β 2 -- Insmed
 TGF- β 3 -- OSI
 Thalassemia gene therapy -- Crucell
 TheraCIM-h-R3 -- Center of Molecular
 Immunology, YM BioSciences
 Theradigm-HBV -- Epimmune
 Theradigm-HPV -- Epimmune
 Theradigm-malaria -- Epimmune
 Theradigm-melanoma -- Epimmune
 TheraFab -- Antisoma
 ThGRF 1-29 -- Theratechnologies
 ThGRF 1-44 -- Theratechnologies
 Thrombin receptor activating peptide --
 Abbott
 thrombomodulin -- Iowa, Novocastra
 Thrombopoietin -- Dragon Pharmaceuticals,
 Genentech
 thrombopoietin, Pliva -- Recepton
 Thrombospondin 2 --
 thrombostatin -- Thromgen
 thymalfasin -- SciClone
 thymocartin -- Gedeon Richter
 thymosin Alpha1 -- NIH
 thyroid stimulating hormone -- Genzyme
 tICAM-1 -- Bayer
 Tick anticoagulant peptide -- Merck
 TIF -- Xoma
 Tifacogin -- Chiron, NIS, Pharmacia
 Tissue factor -- Genentech
 Tissue factor pathway inhibitor
 TJN-135 -- Tsumura
 TM 27 -- Avant
 TM 29 -- Avant
 TMC-151 -- Tanabe Seiyaku
 TNF tumour necrosis factor -- Asahi Kasei
 TNF Alpha -- CytImmune
 TNF antibody -- Johnson & Johnson
 TNF binding protein -- Amgen
 TNF degradation product -- Oncotech
 TNF receptor -- Immunex
 TNF receptor 1, soluble -- Amgen
 TNF Tumour necrosis factor-alpha -- Asahi
 Kasei, Genentech, Mochida
 TNF-Alpha inhibitor -- Tripep
 TNFR:Fc gene therapy -- Targeted Genetics
 TNF-SAM2
 Tolerimab -- Innogenetics
 Toxoplasma gondii vaccine --
 GlaxoSmithKline
 TP 9201 -- Telios
 TP10 -- Avant
 TP20 -- Avant
 tPA -- Centocor
 trafermin -- Scios
 TRAIL/Apo2L -- Immunex
 TRAIL-R1 MAb -- Cambridge Antibody
 Technologies
 transferrin-binding proteins -- CAMR
 Transforming growth factor-beta-1 --
 Genentech
 transport protein -- Genesis
 Trastuzumab -- Genentech
 TRH -- Ferring
 Triabin -- Schering AG
 Triconal
 Triflavin
 troponin I -- Boston Life Sciences
 TRP-2^A -- NIH
 trypsin inhibitor -- Mochida

FIG. 28BB

59/498

<p>TSP-1 gene therapy – TT-232 TTS-CD2 -- Active Biotech Tuberculosis vaccine -- Aventis Pasteur, Genesis Tumor Targeted Superantigens -- Active Biotech -- Pharmacia tumour vaccines -- PhotoCure tumour-activated prodrug antibody conjugates -- Millennium/ImmunoGen tumstatin -- ILEX Tuvirumab -- Novartis TV-4710 -- Teva TWEAK receptor -- Immunex TXU-PAP TY-10721 -- TOA Eiyo Type I diabetes vaccine -- Research Corp Typhoid vaccine CVD 908 U 143677 -- Pharmacia U 81749 -- Pharmacia UA 1248 -- Arizona UGIF -- Sheffield UIC 2 UK 101 UK-279276 -- Corvas Intl. urodilatin -- Pharis urofollitrophin -- Serono Urokinase -- Abbott uteroferrin-- Pepgen V 20 -- GLYCODesign V2 vasopressin receptor gene therapy vaccines -- Active Biotech Varicella zoster glycoprotein vaccine -- Research Corporation Technologies Varicella zoster virus vaccine live -- Cantab Pharmaceuticals Vascular endothelial growth factor -- Genentech, University of California</p>	<p>Vascular endothelial growth factors -- R&D Systems vascular targeting agents -- Peregrine vasopermeation enhancement agents -- Peregrine vasostatin -- NIH VCL -- Bio-Tech. General VEGF -- Genentech, Scios VEGF inhibitor -- Chugai VEGF-2 -- Human Genome Sciences VEGF-Trap -- Regeneron viscumin, recombinant -- Madaus Vitaxin Vitraxe -- ISTA Pharmaceuticals West Nile virus vaccine -- Bavarian Nordic WP 652 WT1 vaccine -- Corixa WX-293 -- Wilex BioTech. WX-360 -- Wilex BioTech. WX-UK1 -- Wilex BioTech. XMP-500 -- XOMA XomaZyme-791 -- XOMA XTL 001 -- XTL Biopharmaceuticals XTL 002 -- XTL Biopharmaceuticals yeast delivery system -- GlobelImmune Yersinia pestis vaccine YIGSR-Stealth -- Johnson & Johnson Yissum Project No. D-0460 -- Yissum YM 207 -- Yamanouchi YM 337 -- Protein Design Labs Yttrium-90 labelled biotin Yttrium-90-labeled anti-CEA MAb T84.66 -- ZD 0490 -- AstraZeneca ziconotide -- Elan ZK 157138 -- Berlex Laboratories Zolimomab arixox Zorcell -- Immune Response ZRXL peptides -- Novartis</p>
---	---

FIG. 28CC